

# **Agilent InfiniiVision 4000 X-Series Oscilloscopes**

## **Service Guide**



**Agilent Technologies**

# Notices

© Agilent Technologies, Inc. 2008-2012

No part of this manual may be reproduced in any form or by any means (including electronic storage and retrieval or translation into a foreign language) without prior agreement and written consent from Agilent Technologies, Inc. as governed by United States and international copyright laws.

## Manual Part Number

54709-97013, November 2012

Print History:

54709-97013, November 2012

54709-97011, October 2012

Available in electronic format only

Agilent Technologies, Inc.  
1900 Garden of the Gods Road  
Colorado Springs, CO 80907 USA

**A newer version of this manual may be available at**

[www.agilent.com/find/4000X-Series-manual](http://www.agilent.com/find/4000X-Series-manual)

Microsoft® is a U.S. registered trademark of Microsoft Corporation.

## Warranty

**The material contained in this document is provided “as is,” and is subject to being changed, without notice, in future editions. Further, to the maximum extent permitted by applicable law, Agilent disclaims all warranties, either express or implied, with regard to this manual and any information contained herein, including but not limited to the implied warranties of merchantability and fitness for a particular purpose. Agilent shall not be liable for errors or for incidental or consequential damages in connection with the furnishing, use, or performance of this document or of any information contained herein. Should Agilent and the user have a separate written agreement with warranty terms covering the material in this document that conflict with these terms, the warranty terms in the separate agreement shall control.**

## Technology Licenses

The hardware and/or software described in this document are furnished under a license and may be used or copied only in accordance with the terms of such license.

## Restricted Rights Legend

If software is for use in the performance of a U.S. Government prime contract or subcontract, Software is delivered and licensed as “Commercial computer software” as defined in DFAR 252.227-7014 (June 1995), or as a “commercial item” as defined in FAR 2.101(a) or as “Restricted computer software” as defined in FAR 52.227-19 (June 1987) or any equivalent agency regulation or contract clause. Use, duplication or disclosure of Software is subject to Agilent Technologies’ standard commercial license terms, and non-DOD Departments and Agencies of the U.S. Government will receive no greater than Restricted Rights as

defined in FAR 52.227-19(c)(1-2) (June 1987). U.S. Government users will receive no greater than Limited Rights as defined in FAR 52.227-14 (June 1987) or DFAR 252.227-7015 (b)(2) (November 1995), as applicable in any technical data.

## Safety Notices

### CAUTION

A **CAUTION** notice denotes a hazard. It calls attention to an operating procedure, practice, or the like that, if not correctly performed or adhered to, could result in damage to the product or loss of important data. Do not proceed beyond a **CAUTION** notice until the indicated conditions are fully understood and met.

### WARNING

A **WARNING** notice denotes a hazard. It calls attention to an operating procedure, practice, or the like that, if not correctly performed or adhered to, could result in personal injury or death. Do not proceed beyond a **WARNING** notice until the indicated conditions are fully understood and met.

# In This Service Guide

This book provides the service information for the Agilent 4000 X-Series oscilloscopes. This manual is divided into these chapters:

## **1 Characteristics and Specifications**

This chapter contains a partial list of characteristics and specifications for the Agilent InfiniiVision 4000 X-Series oscilloscopes.

## **2 Testing Performance**

This chapter explains how to verify correct oscilloscope operation and perform tests to ensure that the oscilloscope meets the performance specifications.

## **3 Calibrating and Adjusting**

This chapter explains how to adjust the oscilloscope for optimum operating performance.

## **4 Troubleshooting**

This chapter begins with suggestions for solving general problems that you may encounter with the oscilloscope. Procedures for troubleshooting the oscilloscope follow the problem solving suggestions.

## **5 Replacing Assemblies**

This chapter describes how to remove assemblies from the 4000 X-Series oscilloscope.

## **6 Replaceable Parts**

This chapter describes how to order replaceable assemblies and parts for the Agilent 4000 X-Series oscilloscopes. It includes diagrams and parts lists for hardware that you can order.

## 7 Safety Notices

At the front of the book you will find safety notice descriptions and document warranties.

### Digital Channels

Because all of the oscilloscopes in the Agilent 4000 X-Series have analog channels, the analog channel topics in this book apply to all instruments. Whenever a topic discusses the digital channels, that information applies only to Mixed-Signal Oscilloscope (MSO) models or DSO models that have been upgraded to an MSO.

### Abbreviated instructions for pressing a series of keys

Instructions for pressing a series of keys are written in an abbreviated manner. Instructions for pressing Key1, then pressing Softkey2, then pressing Softkey3 are abbreviated as follows:

Press [**Key1**] → **Softkey2** → **Softkey3**.

The keys may be front panel keys, or softkeys, which are located directly below the oscilloscope display.



# Contents

- [In This Service Guide](#) 3
- [Figures](#) 9
- [Tables](#) 11
- 1 Characteristics and Specifications** 13
  - Power Requirements 14
  - Measurement Category 15
    - Measurement Category 15
    - Measurement Category Definitions 15
    - Transient Withstand Capability 16
  - Environmental Conditions 17
  - Specifications 18
  - Contact us 19
- 2 Testing Performance** 21
  - Overview 22
  - List of Test Equipment 23
    - Conventions 24
  - To construct the test connector (for use with MSO models only) 25
  - To test digital channels (MSO models only) 27
  - To verify digital channel threshold accuracy (MSO models only) 28
  - When to Test 28

What to Test	28
Verifying Test Results	28
To verify DC vertical gain accuracy	33
To verify dual cursor accuracy	39
To verify bandwidth (-3 dB)	45
To verify time base accuracy	51
To verify trigger sensitivity	53
Test Internal Trigger Sensitivity	54
Test External Trigger Sensitivity (all models)	57
Agilent 4000 X-Series Oscilloscopes Performance Test Record	61
<b>3 Calibrating and Adjusting</b>	<b>65</b>
User Calibration	67
To perform User Cal	67
User Cal Status	70
<b>4 Troubleshooting</b>	<b>71</b>
Solving General Problems with the Oscilloscope	72
If there is no display	72
If there is no trace display	73
If the trace display is unusual or unexpected	73
If you cannot see a channel	74
Verifying Basic Operation	75
To power-on the oscilloscope	75
To perform hardware self test	76
To perform front panel self test	76
To verify default setup	77
To perform an Auto Scale on the Probe Comp signal	78
To compensate passive probes	80
Troubleshooting Internal Assemblies	82

Equipment Required for Troubleshooting Internal Assemblies	82
To prepare for internal assembly troubleshooting	83
Flowchart for Internal Assembly Power Troubleshooting	86
System Board Drawings	87
To check the system board power supply test points	89
To check the line filter board AC output	91
To check the power switch	92
To check the power supply DC output	94
To check the display supplies	95
To check the keyboard supplies	96
To check the fan	97
<b>5 Replacing Assemblies</b>	<b>99</b>
To remove the handle	102
To remove the adjustable legs	103
To remove the cabinet	104
To remove the rear deck assembly	105
To replace the battery	108
To remove the acquisition board	109
To remove the touch controller board	112
To remove the front panel knobs	115
To remove the front bezel assembly	116
To remove the display assembly	120
To remove the keyboard and keypad	122
To remove the fan assembly	125
To remove the power supply shield	126
To remove the line filter board	129

To remove the power supply	132
133	

## **6 Replaceable Parts** 135

Ordering Replaceable Parts	136
Listed Parts	136
Unlisted Parts	136
Direct Mail Order System	136
Exchange Assemblies	137
Exploded Views	138
Replaceable Parts List	141

## **7 Safety Notices** 147

Warnings	147
To clean the instrument	148
Safety Symbols	149

<b>Index</b>	<b>151</b>
--------------	------------

## Figures

Figure 1. Constructing the 8-by-2 Connector	26
Figure 2. Setting Up Equipment for Digital Channel Threshold Accuracy Test	30
Figure 3. Setting up Equipment for DC Vertical Gain Accuracy Test	36
Figure 4. Using a Blocking Capacitor to Reduce Noise	38
Figure 5. Setting up Equipment for Dual Cursor Accuracy Test	42
Figure 6. Using a Blocking Capacitor to Reduce Noise	44
Figure 7. Setting Up Equipment for Bandwidth (-3 dB) Verification Test	47
Figure 8. Setting Up Equipment for Internal Trigger Sensitivity Test	55
Figure 9. Setting Up Equipment for 4-Channel External Trigger Sensitivity Test	58
Figure 10. User Calibration cable for 2-channel oscilloscope	68
Figure 11. User Calibration cable for 4-channel oscilloscope	69
Figure 12. Default setup screen	78
Figure 13. Example pulses	80
Figure 14. Setup for troubleshooting internal assemblies	84
Figure 15. System Board Test Points/Connectors - Top Side, Left	87
Figure 16. System Board Test Points/Connectors - Top Side, Right	88
Figure 17. Verify line filter board AC output	91
Figure 18. Verify power switch operation	93
Figure 19. Location of the Fan Connector	98
Figure 20. Removing the handle	102
Figure 21. Removing adjustable legs	103
Figure 22. Removing the cabinet	104
Figure 23. Removing the BNC securing nuts and washers	105


- Figure 24. Removing the rear deck assembly 106
- Figure 25. Separating front and rear deck assemblies 107
- Figure 26. Battery location 108
- Figure 27. Removing the analog channel BNC securing T6 screws 109
- Figure 28. Acquisition board TORX T10 mounting screw locations 110
- Figure 29. Lifting acquisition board off the front deck 111
- Figure 30. Removing the touch controller shield 112
- Figure 31. Lifting the touch controller shield off 113
- Figure 32. Removing the touch controller board 114
- Figure 33. Removing the front panel knobs 115
- Figure 34. Removing the keyboard cable shield 116
- Figure 35. Disconnecting the keyboard cable 117
- Figure 36. Removing the display cable from the guide 118
- Figure 37. Removing the bezel 119
- Figure 38. Removing the display assembly 120
- Figure 39. Disconnecting the backlight power cable 121
- Figure 40. Disconnecting the softkey board cable 122
- Figure 41. Removing the keyboard 123
- Figure 42. Keyboard and keypad removed 123
- Figure 43. Unclip the softkey board from the bezel 124
- Figure 44. Softkey board and keypad removed 124
- Figure 45. Removing fan assembly 125
- Figure 46. Removing the power supply shield 126
- Figure 47. Power supply shield removed 127
- Figure 48. Removing the air duct 128
- Figure 49. Disconnecting the ground wire 129
- Figure 50. Removing the line filter board 130
- Figure 51. Removing power switch extender 131
- Figure 52. Removing the power supply 132
- Figure 53. Exploded View 1 of 2 139
- Figure 54. Exploded View 2 of 2 140

## Tables

Table 1. List of test equipment	23
Table 2. Conventions	24
Table 3. Materials required to construct the test connectors	25
Table 4. Equipment Required to Test Digital Channel Threshold Accuracy	29
Table 5. Threshold Accuracy Voltage Test Settings	31
Table 6. DC Vertical Gain Accuracy Test Limits	33
Table 7. Equipment Required to Verify DC Vertical Gain Accuracy	34
Table 8. Settings Used to Verify DC Vertical Gain Accuracy	35
Table 9. Equipment Required to Verify Dual Cursor Accuracy	40
Table 10. Settings Used to Verify Dual Cursor Accuracy	41
Table 11. Bandwidth (-3 dB) Test Limits	45
Table 12. Equipment Required to Verify Bandwidth (-3 dB)	46
Table 13. Equipment Required to Verify Time Base Accuracy	51
Table 14. Internal Trigger Sensitivity Test Limits	53
Table 15. External Trigger Sensitivity Test Limits, All Models	53
Table 16. Equipment Required to Verify Trigger Sensitivity	54
Table 17. External Trigger Sensitivity Test Settings	57
Table 18. Equipment Required to Troubleshoot the Oscilloscope	82
Table 19. Replaceable Parts	141





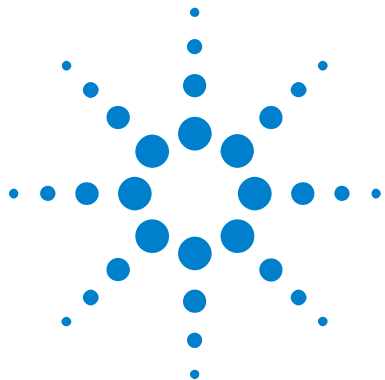

# 1

## Characteristics and Specifications

Power Requirements	14
Measurement Category	15
Environmental Conditions	17
Specifications	18
Contact us	19

This chapter contains a partial list of characteristics and specifications for the Agilent InfiniiVision 4000 X-Series oscilloscopes.

For a full list of Agilent InfiniiVision 4000 X-Series oscilloscopes characteristics and specifications see the data sheets.

The data sheets are available at  
[www.agilent.com/find/4000X-Series](http://www.agilent.com/find/4000X-Series).



## **Power Requirements**

Line voltage, frequency, and power:

- ~Line 100-120 Vac, 50/60/400 Hz
- 100-240 Vac, 50/60 Hz
- 120 W max

## Measurement Category

### Measurement Category

The InfiniiVision 4000 X-Series oscilloscope is intended to be used for measurements in Measurement Category I.

#### **WARNING**

**Use this instrument only for measurements within its specified measurement category.**

---

### Measurement Category Definitions

Measurement category I is for measurements performed on circuits not directly connected to MAINS. Examples are measurements on circuits not derived from MAINS, and specially protected (internal) MAINS derived circuits. In the latter case, transient stresses are variable; for that reason, the transient withstand capability of the equipment is made known to the user.

Measurement category II is for measurements performed on circuits directly connected to the low voltage installation. Examples are measurements on household appliances, portable tools and similar equipment.

Measurement category III is for measurements performed in the building installation. Examples are measurements on distribution boards, circuit-breakers, wiring, including cables, bus-bars, junction boxes, switches, socket-outlets in the fixed installation, and equipment for industrial use and some other equipment, for example, stationary motors with permanent connection to the fixed installation.

Measurement category IV is for measurements performed at the source of the low-voltage installation. Examples are electricity meters and measurements on primary overcurrent protection devices and ripple control units.

## Transient Withstand Capability

### CAUTION

Maximum input voltage for analog inputs

CAT I 300 Vrms, 400 Vpk; transient overvoltage 1.6 kVpk

50  $\Omega$  input: 5 Vrms Input protection is enabled in 50  $\Omega$  mode and the 50  $\Omega$  load will disconnect if greater than 5 Vrms is detected. However the inputs could still be damaged, depending on the time constant of the signal. The 50  $\Omega$  input protection only functions when the oscilloscope is powered on.

With 10073C 10:1 probe: CAT I 500 Vpk, CAT II 400 Vpk

With N2871A, N2872A, N2873A 10:1 probe: CAT I 400 Vpk, transient overvoltage 1.25 kVpk, CAT II 300 Vpk

---

### CAUTION

Maximum input voltage for logic channels:

$\pm 40$  V peak CAT I; transient overvoltage 800 Vpk

---

## Environmental Conditions

Environment	Indoor use only.
Ambient temperature	Operating 0 °C to +55 °C; non-operating –30 °C to +70 °C
Humidity	Operating: 50% to 95% RH at 40 °C for 5 days. Non-operating: 90% RH at 65 °C for 24 hr.
Altitude	Maximum operating altitude: 3,000 m (9,842 ft)
Overvoltage Category	This product is intended to be powered by MAINS that comply to Overvoltage Category II, which is typical of cord-and-plug connected equipment.
Pollution Degree	The InfiniiVision 4000 X-Series oscilloscopes may be operated in environments of Pollution Degree 2 (or Pollution Degree 1).
Pollution Degree Definitions	<p>Pollution Degree 1: No pollution or only dry, non-conductive pollution occurs. The pollution has no influence. Example: A clean room or climate controlled office environment.</p> <p>Pollution Degree 2. Normally only dry non-conductive pollution occurs. Occasionally a temporary conductivity caused by condensation may occur. Example: General indoor environment.</p> <p>Pollution Degree 3: Conductive pollution occurs, or dry, non-conductive pollution occurs which becomes conductive due to condensation which is expected. Example: Sheltered outdoor environment.</p>

## Specifications

Please see the *InfiniiVision 4000 X-Series Oscilloscopes Data Sheet* for complete, up-to-date specifications and characteristics.

To download a copy of the data sheet please visit:  
[www.agilent.com/find/4000X-Series](http://www.agilent.com/find/4000X-Series).

Or go to the Agilent home page at [www.agilent.com](http://www.agilent.com) and search for **4000 X-Series oscilloscopes data sheet**.

To order a data sheet by phone, please contact your local Agilent office. A contact list is provided on the next page. The most up-to-date list is available at:  
[www.agilent.com/find/contactus](http://www.agilent.com/find/contactus)

## Contact us

### Americas

Canada (877) 894-4414  
Latin America 305 269 7500  
United States (800) 829-4444

### Asia Pacific


Australia 1 800 629 485  
China 800 810 0189  
Hong Kong 800 938 693  
India 1 800 112 929  
Japan 81 426 56 7832  
Korea 080 769 0800  
Malaysia 1 800 888 848  
Singapore 1 800 375 8100  
Taiwan 0800 047 866  
Thailand 1 800 226 008

### Europe

Austria 0820 87 44 11  
Belgium 32 (0) 2 404 93 40  
Denmark 45 70 13 15 15  
Finland 358 (0) 10 855 2100  
France 0825 010 700  
Germany 01805 24 6333\*  
\*0.14€/minute  
Ireland 1890 924 204  
Italy 39 02 92 60 8484  
Netherlands 31 (0) 20 547 2111  
Spain 34 (91) 631 3300  
Sweden 0200-88 22 55  
Switzerland (French)  
44 (21) 8113811 (Opt 2)  
Switzerland (German)  
0800 80 53 53 (Opt 1)  
United Kingdom 44 (0) 7004 666666  
Other European countries:  
[www.agilent.com/find/contactus](http://www.agilent.com/find/contactus)





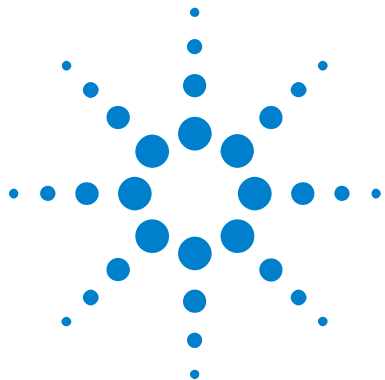

## 2 Testing Performance

Overview	22
List of Test Equipment	23
To construct the test connector (for use with MSO models only)	25
To test digital channels (MSO models only)	27
To verify digital channel threshold accuracy (MSO models only)	28
To verify DC vertical gain accuracy	33
To verify dual cursor accuracy	39
To verify bandwidth (-3 dB)	45
To verify time base accuracy	51
To verify trigger sensitivity	53
Agilent 4000 X-Series Oscilloscopes Performance Test Record	61

This chapter explains how to verify correct oscilloscope operation and perform tests to ensure that the oscilloscope meets the performance specifications.



# Overview

To completely test and troubleshoot MSO models, create and use the test connector accessory as described in this chapter.

- The test connector is only required for oscilloscopes that have the MSO option licensed (enabled).
- The connector is used in the digital channel threshold accuracy test.
- The test connector keeps electrical distortion to a minimum and makes it easy for you to connect the oscilloscope probes to function generators and measurement equipment.

### Let the Equipment Warm Up Before Testing

For accurate test results, let the test equipment and the oscilloscope warm up 30 minutes before testing.

### Verifying Test Results

During the tests, record the readings in the Performance Test Record on [page 61](#). To verify whether a test passes, verify that the reading is within the limits in the Performance Test Record.

**If a performance test fails**

If a performance test fails, first perform the User Cal procedure. Press the following keys to access User Cal: [Utility]→Service→Start User Cal.

## List of Test Equipment

Below is a list of test equipment and accessories required to perform the performance test verification procedures.

**Table 1** List of test equipment

Equipment	Critical Specifications	Recommended Model/ Part Number
Test connector, 8-by-2*	See <a href="#">page 25</a> for instructions on building test connector.	n/a
Digital Multimeter	0.1 mV resolution, 0.005% accuracy	Agilent 34401A
Power Splitter	Outputs differ by 0.15 dB	Agilent 11667B
Oscilloscope Calibrator	DC offset voltage of -5.5 V to 35.5 V, 0.1 V resolution 25 MHz—500 MHz sine wave, 5 ppm	Fluke 5820A
Signal Generator	100 kHz to 1.5 GHz sine waves	Agilent N5181A with 3 GHz option
Power Meter	1.5 GHz $\pm 3\%$ accuracy	Agilent N1914A
Power Sensor	1.5 GHz $\pm 3\%$ accuracy	Agilent E9304A or N8482A
BNC banana cable	BNC (m) to dual banana	Pomona 2BC-BNC-36 or Agilent 11001-66001
BNC cable (qty 3)	BNC - BNC, 48" length	Agilent 10503A <sup>†</sup>
Cable	Type N (m) 609.6 mm (24 in.)	Agilent 11500B
Probe cable*	No substitute	Agilent N6450-60001 (16-channel) or Agilent N6459-60001 (8-channel) <sup>†</sup>
Adapter	BNC(f) to banana(m)	Agilent 1251-2277 <sup>†</sup>
Adapter	BNC Tee (m) (f) (f)	Agilent 1250-0781 <sup>†</sup> or Pomona 3285

\* Required only for testing digital channels of oscilloscopes that have the MSO option.

Most parts and equipment are available at [www.agilent.com](http://www.agilent.com). See respective manufacturer's websites for their equipment.

<sup>†</sup> These parts available at [www.parts.agilent.com](http://www.parts.agilent.com) at the time this manual was published.

**Table 1**    List of test equipment (continued)

Equipment	Critical Specifications	Recommended Model/ Part Number
Adapter	Type N (m) to BNC (m)	Agilent 1250-0082 or Pomona 3288 with Pomona 3533
Blocking capacitor and shorting cap	Note: if a BNC blocking capacitor is not available use an SMA blocking capacitor.	Agilent 11742A + Pomona 4288 + Pomona 5088
Adapter (qty 3)	N(m) to BNC(f)	Agilent 1250-0780
50 Ohm Feedthrough Termination	50Ω BNC (f) to BNC (m)	Agilent 0960-0301

\* Required only for testing digital channels of oscilloscopes that have the MSO option.  
Most parts and equipment are available at [www.agilent.com](http://www.agilent.com). See respective manufacturer's websites for their equipment.  
† These parts available at [www.parts.agilent.com](http://www.parts.agilent.com) at the time this manual was published.

# Conventions

The following conventions will be used when referring to oscilloscope models throughout this chapter.

**Table 2**    Conventions

Models	Referred to as:
MSO-X /DSO-X 4022A, MSO-X /DSO-X 4024A	200 MHz Models
MSO-X /DSO-X 4032A, MSO-X /DSO-X 4034A	350 MHz Models
MSO-X /DSO-X 4052A, MSO-X /DSO-X 4054A	500 MHz Models
MSO-X /DSO-X 4104A	1 GHz Models
MSO-X /DSO-X 4154A	1.5 GHz Models


## To construct the test connector (for use with MSO models only)

Agilent 4000 X-Series oscilloscopes that have digital channels enabled require the test connector described below. Follow the steps to build the test connector.

**Table 3** Materials required to construct the test connectors

Description	Recommended Part	Qty
BNC (f) Connector	Agilent 1250-1032 or Pomona 4578	1
Berg Strip, 8-by-2	3M .100" x .100" Pin Strip Header or similar	1 strip, cut to length (8x2)
Jumper wire		

- 1 Obtain a BNC connector and an 8-by-2 section of Berg strip. A longer strip can be cut to length using wire cutters.
- 2 On one side of the Berg strip, solder a jumper wire to all of the pins (shown in [Figure 1](#) on page 26).
- 3 On the other side of the Berg strip, solder another jumper wire to all of the pins.
- 4 Solder the center of the BNC connector to a center pin on one of the rows on the Berg strip.
- 5 Solder the ground tab of the BNC connector to a center pin on the other row on the Berg strip.

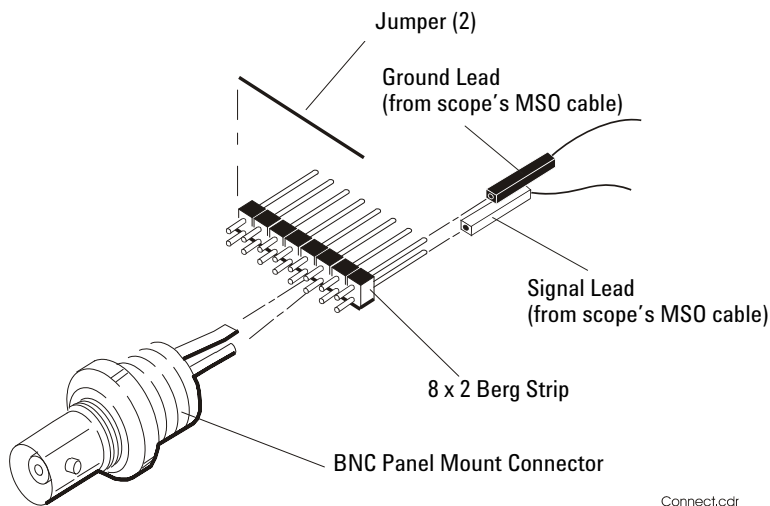

**Figure 1**    Constructing the 8-by-2 Connector

## To test digital channels (MSO models only)

The acquisition system testing provides confidence that the acquisition system is functioning correctly. It does not, however, check a particular specification.

- 1 Disconnect all probes from the circuit under test and from any other input source.
- 2 Using probe leads and grabbers, connect digital channels D0, D1, D2, and D3 to the Probe Comp signal on the center of the front panel.
- 3 Press the [**AutoScale**] key.

If four square waves appear, the acquisition system is functioning correctly.

If the square waves do not appear, go to the “Troubleshooting” chapter. Then return here to finish testing the digital channels.

- 4 Disconnect the digital channels from the calibration point.
- 5 Use steps 2 and 3 to test the following sets of digital channels. After you test one set of digital channels, remove them before connecting the next set.
  - D4, D5, D6, D7
  - D8, D9, D10, D11
  - D12, D13, D14, D15

## To verify digital channel threshold accuracy (MSO models only)

This test verifies the digital channel threshold accuracy specification of the Agilent 4000 X-Series oscilloscopes.

Threshold accuracy test limits:  $\pm(100 \text{ mV} + 3\% \text{ of threshold setting})$

### When to Test

You should perform this test every two years or after 4000 hours of operation, whichever comes first.

### What to Test

Use these instructions to test the threshold settings of digital channels D7-D0. Then, use the same instructions to test digital channels D15-D8.

### Verifying Test Results


After each threshold test, record the voltage reading in the Performance Test Record on [page 61](#). To verify whether a test passes, verify that the voltage reading is within the limits in the Performance Test Record.



**Table 4** Equipment Required to Test Digital Channel Threshold Accuracy

Equipment	Critical Specifications	Recommended Model/Part
Digital Multimeter	0.1 mV resolution, 0.005% accuracy	Agilent 34401A
Oscilloscope Calibrator	DC offset voltage 6.3 V	Fluke 5820A
BNC-Banana Cable		Agilent 11001-66001 or Pomona 2BC-BNC-36
BNC Tee		Agilent 1250-0781 or Pomona 3285
50 $\Omega$ BNC Cable		Agilent 10503A
BNC Test Connector, 8-by-2		User-built (See <a href="#">page 25</a> )
Probe Cable		Agilent N6450-60001 (16-channel) or Agilent N6459-60001 (8-channel)

- 1** Turn on the test equipment and the oscilloscope. Let them warm up for 30 minutes before starting the test.
- 2** Set up the oscilloscope calibrator.
  - a** Set the oscilloscope calibrator to provide a DC offset voltage at the Channel 1 output.
  - b** Use the multimeter to monitor the oscilloscope calibrator DC output voltage.
- 3** Use the 8-by-2 test connector and the BNC cable assembly to connect digital channels D0-D7 to one side of the BNC Tee. Then connect the D0-D7 ground lead to the ground side of the 8-by-2 connector. See [Figure 2](#).

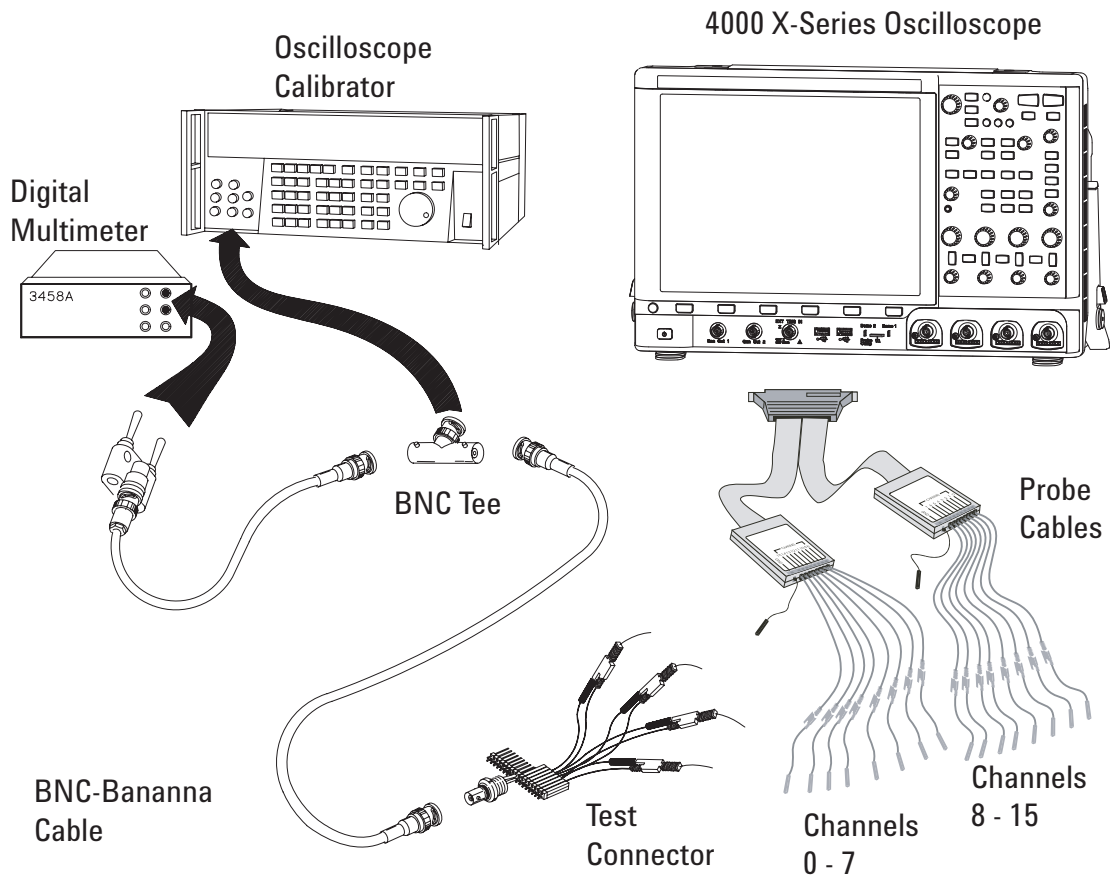

**Figure 2** Setting Up Equipment for Digital Channel Threshold Accuracy Test

- 4 Use a BNC-bananna cable to connect the multimeter to the other side of the BNC Tee.
- 5 Connect the BNC Tee to the Channel 1 output of the calibrator as shown in [Figure 2](#).
- 6 On the oscilloscope, press the **[Digital]** key, then press the **Thresholds** softkey, then press the **D7 - D0** softkey repeatedly until the check mark is next to **User**.

- 7 Press the **User** softkey to the right of the **D7 - D0** softkey, then turn the Entry knob (↻) on the front panel of the oscilloscope to set the threshold test settings as shown in [Table 5](#).

**Table 5** Threshold Accuracy Voltage Test Settings

Threshold voltage setting (in oscilloscope User softkey)	DC offset voltage setting (on oscilloscope calibrator)	Limits
+5.00 V	+5.250 V $\pm$ 1 mV dc	Lower limit = +4.750 V Upper limit = +5.250 V
-5.00 V	-4.750 V $\pm$ 1 mV dc	Lower limit = -5.250 V Upper limit = -4.750 V
0.00 V	+100m V $\pm$ 1 mV dc	Upper limit = +100 mV Lower limit = -100 mV

- 8 Do the following steps for each of the threshold voltage levels shown in [Table 5](#).

- a Set the threshold voltage shown in the **User** softkey using the Entry knob on the oscilloscope.
- b Enter the corresponding DC offset voltage on the oscilloscope calibrator front panel. Then use the multimeter to verify the voltage.

Digital channel activity indicators are displayed on the status line at the top of the oscilloscope display. The activity indicators for D7-D0 should show all of the channels at digital high levels.

- c Use the knob on the oscilloscope calibrator to decrease the offset voltage, in increments of 10 mV, until the activity indicators for digital channels D7-D0 are all at digital low levels. Record the oscilloscope calibrator voltage in the Performance Test Record (see [page 61](#)).
- d Use the knob on the oscilloscope calibrator to increase the offset voltage, in increments of 10 mV, until the activity indicators for digital channels D7-D0 are all at

digital high levels. Record the oscilloscope calibrator voltage in the Performance Test Record (see [page 61](#)).

Before proceeding to the next step, make sure that you have recorded the oscilloscope calibrator voltage levels for each of the threshold settings shown in [Table 5](#).

- 9 When testing 4000 X-Series MSOs, use the 8-by-2 test connector to connect digital channels D15-D8 to the output of the oscilloscope calibrator. Then connect the D15-D8 ground lead to the ground side of the 8-by-2 connector.
- 10 Repeat this procedure (steps 6 through 8) for digital channels D15-D8 to verify threshold accuracy and record the threshold levels in the Performance Test Record (see [page 61](#)). Be sure to set the thresholds with the **User** softkey for the appropriate set of channels.

## To verify DC vertical gain accuracy

This test verifies the accuracy of the analog channel DC vertical gain for each channel.

In this test, you will measure the dc voltage output of an oscilloscope calibrator using the oscilloscope's

**Average - Full Screen** voltage measurement and compare the results with the multimeter reading.

**Table 6** DC Vertical Gain Accuracy Test Limits

Models	Test Limits	Notes
4000 X-Series	±2.0% of full scale	<ul style="list-style-type: none"><li>• Full scale is defined as 32 mV on the 2 mV/div range and the 1 mV/div range.</li><li>• Full scale on all other ranges is defined as 8 divisions times the V/div setting.</li></ul>

**Table 7**    Equipment Required to Verify DC Vertical Gain Accuracy

Equipment	Critical Specifications	Recommended Model/Part
Oscilloscope Calibrator	14 mV to 35 Vdc, 0.1 V resolution	Fluke 5820A
Digital multimeter	Better than 0.01% accuracy	Agilent 34401A
Cable	BNC, Qty 2	Agilent 10503A
Shorting cap	BNC	Agilent 1250-0774
Adapter	BNC (f) to banana (m)	Agilent 1251-2277
Adapter	BNC tee (m) (f) (f)	Agilent 1250-0781 or Pomona 3285
Blocking capacitor		Agilent 11742A + Pomona 4288 + Pomona 5088

- 1 Press **[Save/Recall] > Default/Erase > Factory Default** to recall the factory default setup.
- 2 Set up the oscilloscope.
  - a Adjust the horizontal scale to **10.00 ms/div**.
  - b Set the Volts/Div setting to the value in the first line in [Table 8](#) (depending on the oscilloscope model).
  - c Adjust the channel's vertical position knob to place the baseline (reference level) at 0.5 major division from the bottom of the display.

**Table 8** Settings Used to Verify DC Vertical Gain Accuracy

Volts/Div Setting	Oscilloscope Calibrator Setting	Test Limits		
5 V/Div	35 V	34.2 V	to	35.8 V
2 V/Div	14 V	13.68 V	to	14.32 V
1 V/Div	7 V	6.84 V	to	7.16 V
500 mV/Div	3.5 V	3.42 V	to	3.58 V
200 mV/Div	1.4 V	1.368 V	to	1.432 V
100 mV/Div	700 mV	684 mV	to	716 mV
50 mV/Div	350 mV	342 mV	to	358 mV
20 mV/Div	140 mV	136.8 mV	to	143.2 mV
10 mV/Div	70 mV	68.4 mV	to	71.6 mV
5 mV/Div <sup>1</sup>	35 mV	34.2 mV	to	35.8 mV
2 mV/Div <sup>1,2</sup>	14 mV	13.36 mV	to	14.64 mV
1 mV/Div <sup>1,2</sup>	7 mV	6.36 mV	to	7.64 mV

<sup>1</sup> A blocking capacitor is required at this range to reduce noise. See [“Use a Blocking Capacitor to Reduce Noise”](#) on page 38.

<sup>2</sup> Full scale is defined as 32 mV on the 2 mV/div range and the 1 mV/div range. Full scale on all other ranges is defined as 8 divisions times the V/div setting.

**d** Press the **[Acquire]** key.


**e** Then press the **Acq Mode** softkey and select **Averaging**.

**f** Then press the **#Avgs** softkey and set it to 64.

Wait a few seconds for the measurement to settle.

**3** Add a measurement for the average voltage:

**a** Press the **[Meas]** key.

**b** Press **Source**; then, turn the Entry knob (labeled  on the front panel) to select the channel you are testing.





- 6** Adjust the output so that the multimeter reading displays the first Volts/div calibrator setting value in [Table 8](#) (depending on the oscilloscope model).
- 7** Disconnect the multimeter.
- 8** Wait until the measurement settles.
- 9** Read the “current” average voltage value again as V2.
- 10** Calculate the difference  $V2 - V1$ .


The difference in average voltage readings should be within the test limits of [Table 8](#) (depending on the oscilloscope model).

If a result is not within the test limits, go to the “Troubleshooting” chapter. Then return here.

- 11** Disconnect the oscilloscope calibrator from the oscilloscope.
- 12** Repeat this procedure to check the DC vertical gain accuracy with the remaining Volts/div setting values in [Table 8](#) (depending on the oscilloscope model).
- 13** Finally, repeat this procedure for the remaining channels to be tested.

### Use a Blocking Capacitor to Reduce Noise

On the more sensitive ranges, such as 1 mV/div, 2 mV/div, and 5 mV/div, noise may be a factor. To eliminate the noise, add a BNC Tee, blocking capacitor, and shorting cap at the oscilloscope channel input to shunt the noise to ground. See [Figure 4](#). If a BNC capacitor is not available, use an SMA blocking capacitor, adapter, and cap. See [“Blocking capacitor and shorting cap”](#) in the equipment list on [page 24](#) for details.

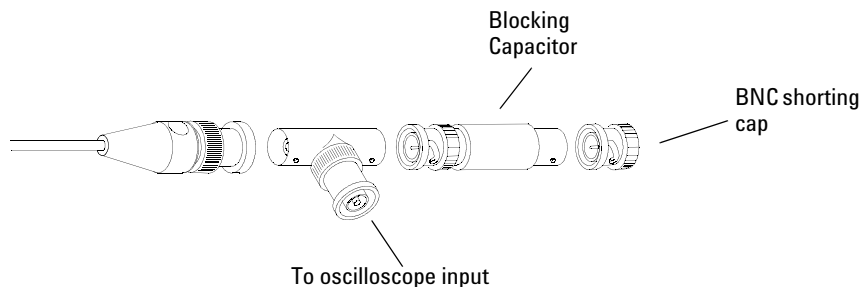

**Figure 4** Using a Blocking Capacitor to Reduce Noise

## To verify dual cursor accuracy

This test verifies the dual cursor accuracy for each analog channel.

This test is similar to the test for verifying the DC vertical gain, except you will measure the dc voltage output of an oscilloscope calibrator using dual cursors on the oscilloscope and compare the results with the multimeter reading.

Dual cursor accuracy test limits:  $\pm$ [DC vertical gain accuracy + 0.42% full scale]

For the DC vertical gain accuracy test limits, see [Table 6](#) on page 33.

**Table 9**    Equipment Required to Verify Dual Cursor Accuracy

<b>Equipment</b>	<b>Critical Specifications</b>	<b>Recommended Model/Part</b>
Oscilloscope Calibrator	14 mV to 35 Vdc, 0.1 V resolution	Fluke 5820A
Digital multimeter	Better than 0.01% accuracy	Agilent 34401A
Cable	BNC, Qty 2	Agilent 10503A
Shorting cap	BNC	Agilent 1250-0774
Adapter	BNC (f) to banana (m)	Agilent 1251-2277
Adapter	BNC tee (m) (f) (f)	Agilent 1250-0781 or Pomona 3285
Blocking capacitor		Agilent 11742A + Pomona 4288 + Pomona 5088

- 1** Press **[Save/Recall] > Default/Erase > Factory Default** to recall the factory default setup.
- 2** Set up the oscilloscope.
  - a** Set the Volts/Div setting to the value in the first line in [Table 10](#) (depending on the oscilloscope model).
  - b** Adjust the channel 1 position knob to place the baseline at 0.5 major division from the bottom of the display.

**Table 10** Settings Used to Verify Dual Cursor Accuracy


Volts/Div Setting	Oscilloscope Calibrator Setting	Test Limits		
5 V/Div	35 V	34.032 V	to	35.968 V
2 V/Div	14 V	13.6128 V	to	14.3872 V
1 V/Div	7 V	6.8064 V	to	7.1936 V
500 mV/Div	3.5 V	3.4032 V	to	3.5968 V
200 mV/Div	1.4 V	1.36128 V	to	1.43872 V
100 mV/Div	700 mV	680.64 mV	to	719.36 mV
50 mV/Div	350 mV	340.32 mV	to	359.68 mV
20 mV/Div	140 mV	136.128 mV	to	143.87 mV
10 mV/Div	70 mV	68.064 mV	to	71.94 mV
5 mV/Div <sup>1</sup>	35 mV	34.032 mV	to	35.97 mV
2 mV/Div <sup>1,2</sup>	14 mV	13.226 mV	to	14.77 mV
1 mV/Div <sup>1,2</sup>	7 mV	6.226 mV	to	7.77 mV

<sup>1</sup> A blocking capacitor is required at this range to reduce noise. See [“Use a Blocking Capacitor to Reduce Noise”](#) on page 44.


<sup>2</sup> Full scale is defined as 32 mV on the 2 mV/div range and the 1 mV/div range. Full scale on all other ranges is defined as 8 divisions times the V/div setting.

- c** Press the **[Acquire]** key.
- d** Then press the **Acq Mode** softkey and select **Averaging**.
- e** Then press the **#Avgs** softkey and set it to 64.

Wait a few seconds for the measurement to settle.

- 3** Press the **[Cursors]** key, set the **Mode** softkey to **Normal**, then press the **XY** softkey and select **Y**. Press the **Y1** softkey, then use the Entry knob (labeled  on the front panel) to set the Y1 cursor on the baseline of the signal.

- 4 Use the BNC tee and cables to connect the oscilloscope calibrator /power supply to both the oscilloscope and the multimeter (see [Figure 5](#)).

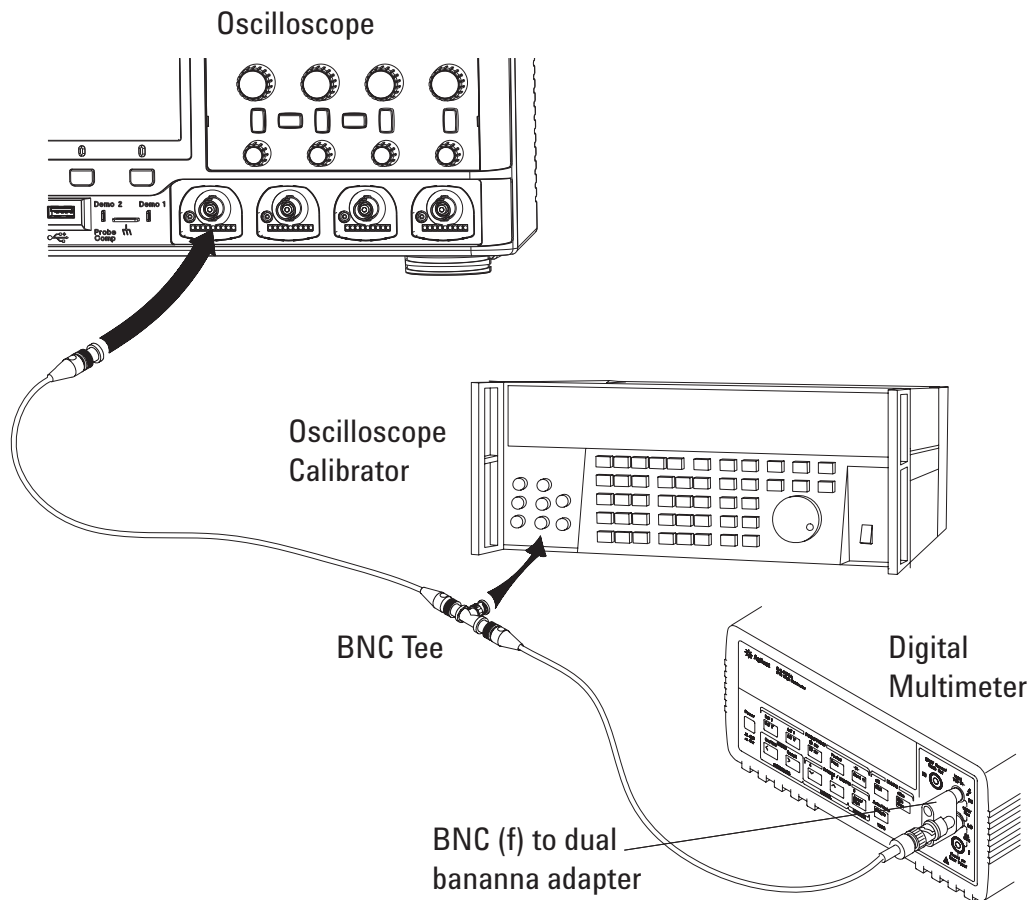

**Figure 5** Setting up Equipment for Dual Cursor Accuracy Test

- 5 Adjust the output so that the multimeter reading displays the first Volts/div calibrator setting value in [Table 10](#).
- 6 Disconnect the multimeter.
- 7 Wait until the measurement settles.

- 8 Press the **Y2** softkey, then position the Y2 cursor to the center of the voltage trace using the Entry knob.


The  $\Delta Y$  value on the lower line of the display should be within the test limits of [Table 10](#).

If a result is not within the test limits, go to the “Troubleshooting” chapter. Then return here.

- 9 Disconnect the oscilloscope calibrator from the oscilloscope.
- 10 Repeat this procedure to check the dual cursor accuracy with the remaining Volts/div setting values in [Table 10](#).
- 11 Finally, repeat this procedure for the remaining channels to be tested.

### Use a Blocking Capacitor to Reduce Noise

On the more sensitive ranges, such as 1 mV/div, 2 mV/div, and 5 mV/div, noise may be a factor. To eliminate the noise, add a BNC Tee, blocking capacitor, and shorting cap at the oscilloscope channel input to shunt the noise to ground. See [Figure 6](#). If a BNC capacitor is not available, use an SMA blocking capacitor, adapter, and cap. See [“Blocking capacitor and shorting cap”](#) in the equipment list on [page 24](#) for details.

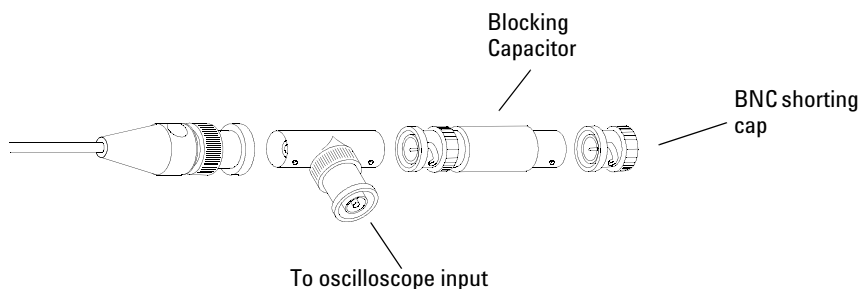

**Figure 6** Using a Blocking Capacitor to Reduce Noise



## To verify bandwidth (-3 dB)

This test checks the bandwidth (-3 dB) of the oscilloscope. In this test you will use a signal generator and a power meter.


**Table 11** Bandwidth (-3 dB) Test Limits

Models	Test Limits
1.5 GHz Models	All channels (-3 dB), dc to 1.5 GHz
1 GHz Models	All channels (-3 dB), dc to 1 GHz
500 MHz Models	All channels (-3 dB), dc to 500 MHz
350 MHz Models	All channels (-3 dB), dc to 350 MHz
200 MHz Models	All channels (-3 dB), dc to 200 MHz

**Table 12**    Equipment Required to Verify Bandwidth (-3 dB)

Equipment	Critical Specifications	Recommended Model/Part
Signal Generator	100 kHz - 1.5 GHz at 200 mVrms	Agilent N5181A with 3 GHz option
Power Meter	1 MHz - 1.5 GHz $\pm 3\%$ accuracy	Agilent N1914A
Power Sensor	1 MHz - 1.5 GHz $\pm 3\%$ accuracy	Agilent E9304A or N8482A
Power Splitter	outputs differ by $< 0.15$ dB	Agilent 11667A
Cable	Type N (m) 24 inch	Agilent 11500B
Adapter	Type N (m) to BNC (m)	Agilent 1250-0082 or Pomona 3288 with Pomona 3533

- 1** Connect the equipment (see [Figure 7](#)).
  - a** Use the N cable to connect the signal generator to the input of the power splitter input.
  - b** Connect the power sensor to one output of the power splitter.
  - c** Use an N-to-BNC adapter to connect the other splitter output to the channel 1 input.

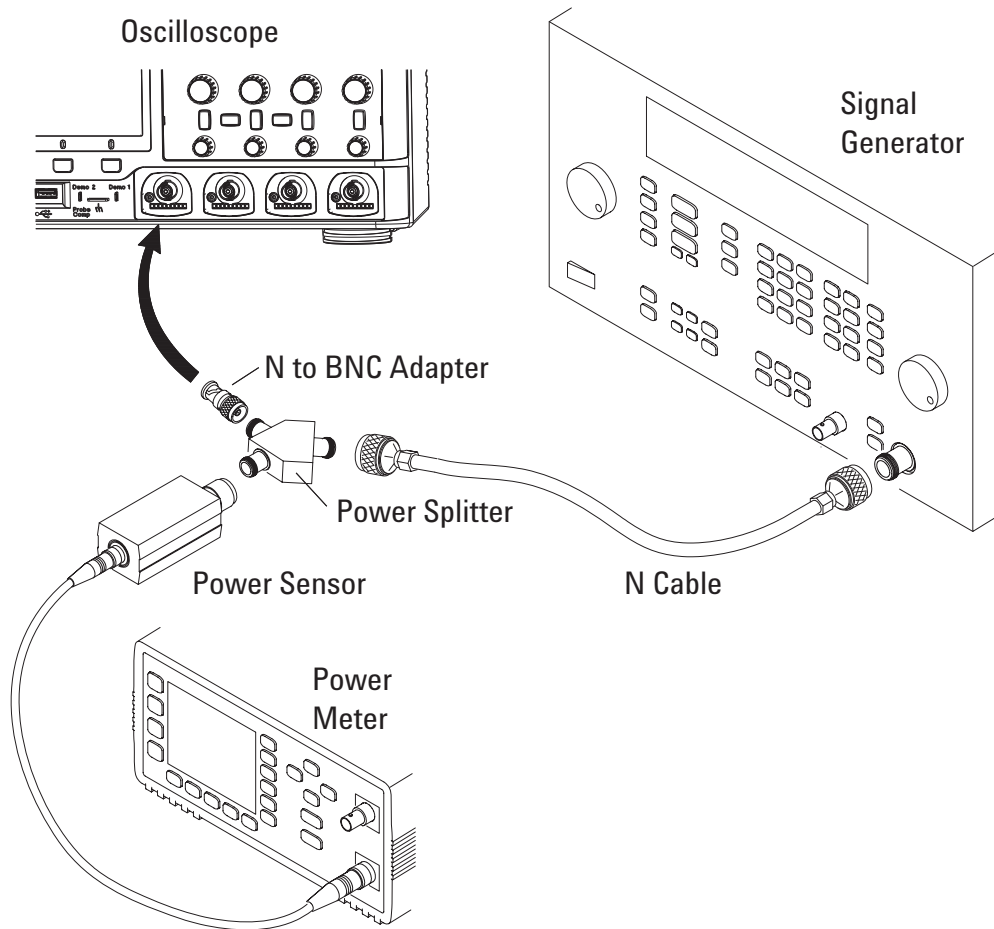

**Figure 7** Setting Up Equipment for Bandwidth (-3 dB) Verification Test

**2** Set up the power meter.

Set the power meter to display measurements in units of watts.

**3** Set up the oscilloscope.

- a** Press the **[Default Setup]** key.
- b** Set channel 1 **Coupling** to **DC**.
- c** Set channel 1 **Imped** to **50 Ohm**.

- d** Set the time base to 500 ns/div.
  - e** Set the Volts/Div for channel 1 to 200 mV/div.
  - f** Press the **[Acquire]** key, then press the **Averaging** softkey.
  - g** Turn the Entry knob to set **# Avgs** to 8 averages.
- 4** Set the signal generator for 1 MHz and six divisions of amplitude.

The signal on the oscilloscope screen should be about five cycles at six divisions amplitude.

- 5** Set up the Amplitude measurement
- a** Press the **[Meas]** key.
  - b** Press the **Clear Meas** softkey and then the **Clear All** softkey.
  - c** Press the **Type:** softkey and use the Entry knob to select **AC RMS - Full Screen (Std Deviation)** within the select menu.
  - d** Press the **Add Measurement** softkey.
- 6** Note the oscilloscope AC RMS - FS(1) reading at the bottom of the screen. (This is the RMS value with any dc offset removed.)
- 7** Set the power meter Cal Factor % to the 1 MHz value on the calibration chart on the power sensor.
- 8** Note the reading on the power meter and covert to Vrms using the expression:

$$V_{in_{1MHz}} = \sqrt{P_{meas_{1MHz}} \times 50\Omega}$$

For example, if the power meter reading is 892 uW, then  
 $V_{in_{1MHz}} = (892 \times 10^{-6} \times 50\Omega)^{1/2} = 211.2 \text{ mV}_{\text{rms}}$ .

- 9** Change the signal generator output frequency according to the maximum frequency for the oscilloscope using the following:
- 1.5 GHz Models: 1.5 GHz
  - 1 GHz Models: 1 GHz
  - 500 MHz Models: 500 MHz
  - 350 MHz Models: 350 MHz
  - 200 MHz Models: 200 MHz
- 10** Referencing the frequency from step 9, set the power meter Cal Factor % to the frequency value on the calibration chart on the power sensor.
- 11** Set the oscilloscope sweep speed according to the following:
- 1.5 GHz Models: 500 ps/div
  - 1 GHz Models: 500 ps/div
  - 500 MHz Models: 1 ns/div
  - 350 MHz Models: 2 ns/div
  - 200 MHz Models: 2 ns/div
- 12** Note the oscilloscope AC RMS - FS(1) reading at the bottom of the screen.
- 13** Note the reading on the power meter and covert to Vrms using the expression:

$$Vin_{maxfreq} = \sqrt{Pmeas_{maxfreq} \times 50\Omega}$$

- 14** Calculate the response using the expression:

$$\text{response(dB)} = 20 \log_{10} \left[ \frac{Vout_{\max \text{ freq}} / Vin_{\max \text{ freq}}}{Vout_{1\text{MHz}} / Vin_{1\text{MHz}}} \right]$$

**Example**

If:

$$Pmeas_{1\text{MHz}} = 892 \text{ uW}$$

$$\text{AC RMS - FS(n)}_{1\text{MHz}} = 210.4 \text{ mV}$$

$$Pmeas_{\max \text{ freq}} = 687 \text{ uW}$$

$$\text{AC RMS - FS(n)}_{\max \text{ freq}} = 161.6 \text{ mV}$$

Then after converting the values from the power meter to  $V_{rms}$ :

$$\text{response(dB)} = 20 \log_{10} \left[ \frac{161.6 \text{ mV} / 185.3 \text{ mV}}{210.4 \text{ mV} / 211.2 \text{ mV}} \right] = -1.16 \text{ dB}$$

- 15** The result from step 14 should be within -3.0 dB. Record the result in the Performance Test Record (see [page 61](#)).
- 16** Move the power splitter from the channel 1 to the channel 2 input.
- 17** Turn off the current channel and turn on the next channel using the channel keys.
- 18** Repeat steps 3 through 17 for the remaining channels, setting the parameters of the channel being tested where appropriate.

## To verify time base accuracy

This test verifies the accuracy of the time base. In this test you will measure the absolute error of the time base oscillator and compare the results to the specification.

**Table 13** Equipment Required to Verify Time Base Accuracy

Equipment	Critical Specifications	Recommended Model/Part
Signal Generator	100 kHz - 1.5 GHz, 0.01 Hz frequency resolution, jitter: < 2ps	Agilent N5181A with 3 GHz option
Cable	BNC, 3 feet	Agilent 10503A

- 1 Set up the signal generator.
  - a Set the output to 10 MHz, approximately 1 V<sub>pp</sub> sine wave.
- 2 Connect the output of the signal generator to oscilloscope channel 1 using the BNC cable.
- 3 Set up the oscilloscope:
  - a Press [AutoScale].
  - b Set the oscilloscope Channel 1 vertical sensitivity to 200 mv/div.
  - c Set the oscilloscope horizontal sweep speed control to 5 ns/div.
  - d Adjust the intensity to get a sharp, clear trace.
  - e Adjust the oscilloscope's trigger level so that the rising edge of the waveform at the center of the screen is located where the center horizontal and vertical grid lines cross (center screen).
  - f Ensure the horizontal position control is set to 0.0 seconds.

- 4 Make the measurement.
  - a Set oscilloscope horizontal sweep speed control to 1 ms/div.
  - b Set horizontal position control to +1 ms (rotate control CCW).
  - c Set the oscilloscope horizontal sweep speed control to 5 ns/div.
  - d Record the number of nanoseconds from where the rising edge crosses the center horizontal grid line to the center vertical grid line. The number of nanoseconds is equivalent to the time base error in ppm.

Time base accuracy limit:  $\pm 10$  ppm

- e Record the result and compare it to the limits in the Performance Test Record (see [page 61](#)).



## To verify trigger sensitivity

This test verifies the trigger sensitivity. In this test, you will apply a sine wave to the oscilloscope at the upper bandwidth limit. You will then decrease the amplitude of the signal to the specified levels, and check to see if the oscilloscope is still triggered.

**Table 14** Internal Trigger Sensitivity Test Limits

BW Models	V/div	Frequency	Sensitivity
200 MHz, 350 MHz, 500 MHz, and 1 GHz	< 10 mV/div	All	greater of 1 div or 5 mV <sub>pp</sub>
	>= 10 mV/div	All	0.6 div
1.5 GHz	< 10 mV/div	DC to 1 GHz	greater of 1 div or 5 mV <sub>pp</sub>
		1 GHz to 1.5 GHz	greater of 1.5 div or 5 mV <sub>pp</sub>
	>= 10 mV/div	DC to 1 GHz	0.6 div
		1 GHz to 1.5 GHz	1.0 div

**Table 15** External Trigger Sensitivity Test Limits, All Models


Input Range	Frequency	Sensitivity
1.6 V	DC to 100 MHz	40 mV <sub>pp</sub>
	100 MHz to 200 MHz	70 mV <sub>pp</sub>
8 V	DC to 100 MHz	200 mV <sub>pp</sub>
	100 MHz to 200 MHz	350 mV <sub>pp</sub>

**Table 16**    Equipment Required to Verify Trigger Sensitivity

Equipment	Critical Specifications	Recommended Model/Part
Signal Generator	100 kHz to 1.5 GHz sine waves	Agilent N5181A with 3 GHz option
Power splitter	Outputs differ < 0.15 dB	Agilent 11667A
Power Meter	1.5 GHz $\pm 3\%$ accuracy	Agilent N1914A
Power Sensor	1.5 GHz $\pm 3\%$ accuracy	Agilent E9304A or N8482A
Cable	BNC, Qty 3	Agilent 10503A
Adapter	N (m) to BNC (f), Qty 3	Agilent 1250-0780
Feedthrough	50 $\Omega$ BNC (f) to BNC (m)	Agilent 0960-0301

## Test Internal Trigger Sensitivity

- 1 On the oscilloscope, press the **[Default Setup]** key.
- 2 Press the **[Mode/Coupling]** key; then, press the **Mode** softkey to select **Normal**.
- 3 Connect the equipment (see [Figure 8](#)).
  - a Connect the signal generator output to the oscilloscope channel 1 input.

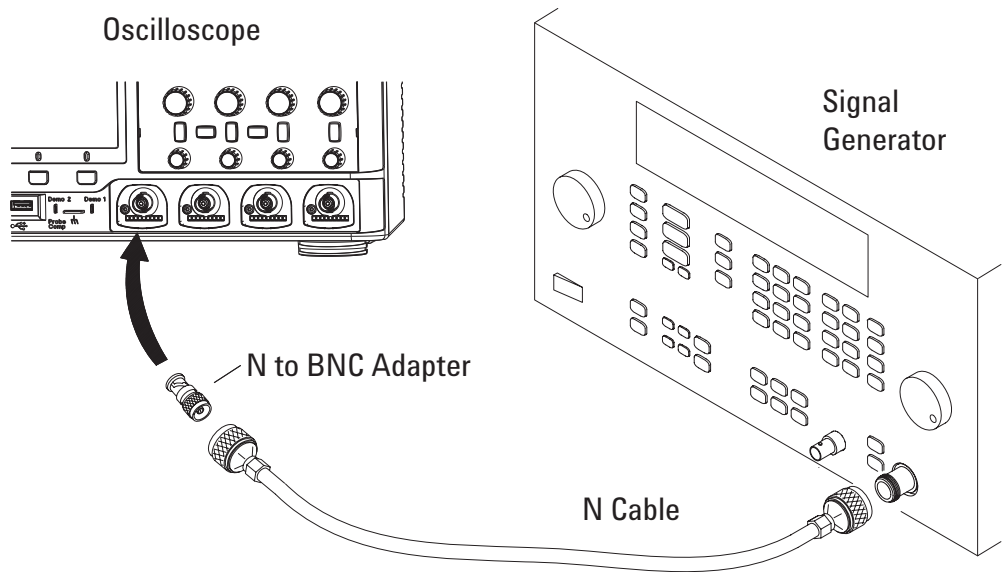

**Figure 8** Setting Up Equipment for Internal Trigger Sensitivity Test

- b** Set channel 1 **Imped** to **50 Ohm**.
- 4** To verify the trigger sensitivity at the oscilloscope's maximum bandwidth, set the output frequency of the signal generator to the maximum bandwidth of the oscilloscope:
  - 1.5 GHz models: 1.5 GHz and 1 GHz.
  - 1 GHz models: 1 GHz.
  - 500 MHz models: 500 MHz.
  - 350 MHz models: 350 MHz.
  - 200 MHz models: 200 MHz.
- 5** Perform these steps to test at the 5 mV/div setting:
  - a** Set the signal generator amplitude to about 10 mV<sub>pp</sub>.
  - b** Press the **[AutoScale]** key.
  - c** Set the time base to 10 ns/div.
  - d** Set channel 1 to 5 mV/div.

- e** Decrease the amplitude from the signal generator until 1 vertical division of the signal (about 5 mV<sub>pp</sub>) is displayed.

The trigger is stable when the displayed waveform is stable. If the trigger is not stable, try adjusting the trigger level. If adjusting the trigger level makes the trigger stable, the test still passes. If adjusting the trigger does not help, see the “Troubleshooting” chapter. Then return here.

- f** Record the result as Pass or Fail in the Performance Test Record (see [page 61](#)).
  - g** Repeat this step for the remaining oscilloscope channels.
- 6** Perform these steps to test at the 10 mV/div setting:
- a** Set the signal generator amplitude to about 20 mV<sub>pp</sub>.
  - b** Press the **[AutoScale]** key.
  - c** Set the time base to 10 ns/div.
  - d** Set channel 1 to 10 mV/div.
  - e** Decrease the amplitude from the signal generator until 0.6 vertical divisions of the signal (about 6 mV<sub>pp</sub>) is displayed.

The trigger is stable when the displayed waveform is stable. If the trigger is not stable, try adjusting the trigger level. If adjusting the trigger level makes the trigger stable, the test still passes. If adjusting the trigger does not help, see the “Troubleshooting” chapter. Then return here.

- f** Record the result as Pass or Fail in the Performance Test Record (see [page 61](#)).
- g** Repeat this step for the remaining oscilloscope channels.

## Test External Trigger Sensitivity (all models)


This test applies to all models.

Verify the external trigger sensitivity at these settings:

**Table 17** External Trigger Sensitivity Test Settings

Input Range	100 MHz	200 Mhz
1.6 V	40 mV <sub>pp</sub>	70 mV <sub>pp</sub>
8 V	200 mV <sub>pp</sub>	350 mV <sub>pp</sub>

- 1 Connect the equipment (see [Figure 9](#)).
  - a Use the N cable to connect the signal generator to the power splitter input.
  - b Connect one output of the power splitter to the Aux Trig input through a 50Ω feedthrough termination.
  - c Connect the power sensor to the other output of the power splitter.

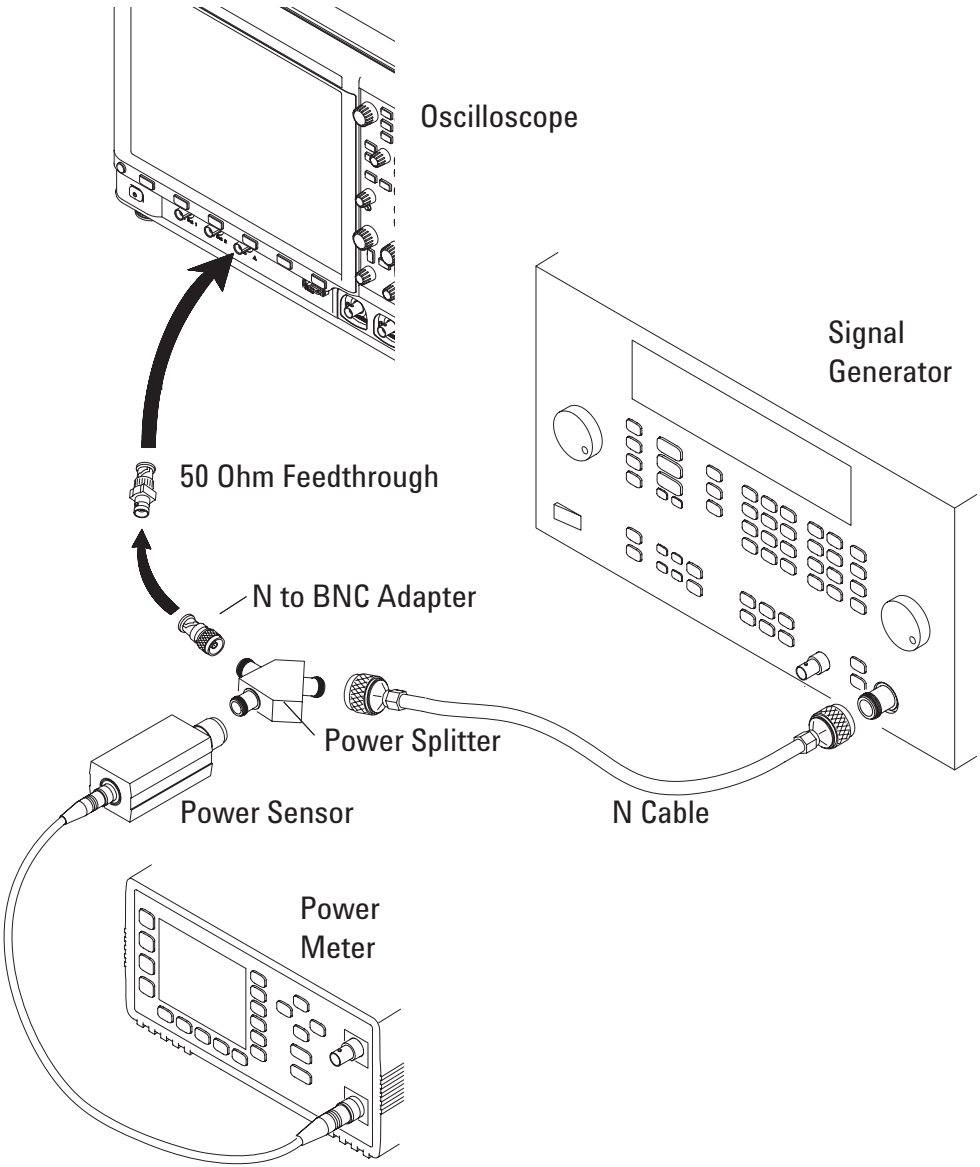

**Figure 9**    Setting Up Equipment for 4-Channel External Trigger Sensitivity Test

- 2 Set up the oscilloscope.
  - a Press the **[Default Setup]** key.
  - b Press the **[Mode/Coupling]** key; then, press the **Mode** softkey to select **Normal**.
- 3 Change the signal generator output frequency to 100 MHz or 200 MHz.
- 4 Set the power meter Cal Factor % to the appropriate value (100 MHz or 200 MHz) on the calibration chart on the power sensor. If necessary, do a linear interpolation if a 100 MHz or 200 MHz factor is not included in the power meter's calibration chart.
- 5 Adjust the signal generator output for reading on the power meter of:

Signal Generator Frequency	Calculation	Power Meter Reading
100 MHz	$40 \text{ mV}_{pp} = 14.14 \text{ mV rms,}$ $\text{Power} = V_{in}^2/50\Omega = 14.14 \text{ mV}^2/50\Omega$	4 $\mu\text{W}$
200 MHz	$70 \text{ mV}_{pp} = 24.75 \text{ mV rms,}$ $\text{Power} = V_{in}^2/50\Omega = 24.75 \text{ mV}^2/50\Omega$	12.25 $\mu\text{W}$
100 MHz	$200 \text{ mV}_{pp} = 70.71 \text{ mV rms,}$ $\text{Power} = V_{in}^2/50\Omega = 70.71 \text{ mV}^2/50\Omega$	100 $\mu\text{W}$
200 MHz	$350 \text{ mV}_{pp} = 123.74 \text{ mV rms,}$ $\text{Power} = V_{in}^2/50\Omega = 123.74 \text{ mV}^2/50\Omega$	306 $\mu\text{W}$

- 6 Press the **[Trigger]** key, then press the **Source** softkey to set the trigger source to **External**.
- 7 Check for stable triggering and adjust the trigger level if necessary. Triggering is indicated by the **Trig'd** indicator at the top of the display. When **Trig'd?** is displayed, the oscilloscope is not triggered. When **Trig'd** is displayed (no question mark), the oscilloscope is triggered.

- 8 Record the results as Pass or Fail in the Performance Test Record (see [page 61](#)).

If the test fails, see the “Troubleshooting” chapter. Then return here.



# Agilent 4000 X-Series Oscilloscopes Performance Test Record

Serial No. _____			Test by _____	
Test Interval _____			Work Order No. _____	
Recommended Next Testing _____			Temperature _____	

Threshold	Specification	Limits	Ch D7-D0	Ch D15-D8
<b>Accuracy Test</b>	5 V - 250 mV	4.750 V	_____	_____
<b>(100 mV + 3% of threshold setting)</b>	5 V + 250 mV	5.250 V	_____	_____
	-5 V - 250 mV	-5.250 V	_____	_____
	-5 V + 250 mV	-4.750 V	_____	_____
	0 V - 100 mV	-100 mV	_____	_____
	0 V + 100 mV	100 mV	_____	_____

DC Vertical Gain Accuracy			Channel 1	Channel 2	Channel 3*	Channel 4*
Range	Power Supply Setting	Test Limits				
5 V/Div	35 V	34.2 V to 35.8 V	_____	_____	_____	_____
2 V/Div	14 V	13.68 V to 14.32 V	_____	_____	_____	_____
1 V/Div	7 V	6.84 V to 7.16 V	_____	_____	_____	_____
500 mV/Div	3.5 V	3.42 V to 3.58 V	_____	_____	_____	_____
200 mV/Div	1.4 V	1.368 V to 1.432 V	_____	_____	_____	_____
100 mV/Div	700 mV	684 mV to 716 mV	_____	_____	_____	_____
50 mV/Div	350 mV	342 mV to 358 mV	_____	_____	_____	_____
20 mV/Div	140 mV	136.8 mV to 143.2 mV	_____	_____	_____	_____
10 mV/Div	70 mV	68.4 mV to 71.6 mV	_____	_____	_____	_____
5 mV/Div	35 mV	34.2 mV to 35.8 mV	_____	_____	_____	_____
2 mV/Div	14 mV	13.36 mV to 14.64 mV	_____	_____	_____	_____
1 mV/Div	7 mV	6.36 mV to 7.64 mV	_____	_____	_____	_____

Continued on next page.

## 2 Testing Performance

Dual Cursor Accuracy				Channel 1	Channel 2	Channel 3*	Channel 4*
Range	Power Supply Setting	Test Limits					
5 V/Div	35 V	34.0 V to 36.0 V					
2 V/Div	14 V	13.6 V to 14.4 V					
1 V/Div	7 V	6.8 V to 7.2 V					
500 mV/Div	3.5 V	3.4 V to 3.6 V					
200 mV/Div	1.4 V	1.36 V to 1.44 V					
100 mV/Div	700 mV	680 mV to 720 mV					
50 mV/Div	350 mV	340 mV to 360 mV					
20 mV/Div	140 mV	136 mV to 144 mV					
10 mV/Div	70 mV	68 mV to 72 mV					
5 mV/Div	35 mV	34 mV to 36 mV					
2 mV/Div	14 mV	13.2 mV to 14.8 mV					
1 mV/Div	7 mV	6.2 mV to 7.8 mV					
Bandwidth (-3 dB)	Model	Test Limits		Channel 1	Channel 2	Channel 3*	Channel 4*
	415x	-3 dB at 1.5 GHz					
	410x	-3 dB at 1 GHz					
	405x	-3 dB at 500 MHz					
	403x	-3 dB at 350 MHz					
	402x	-3 dB at 200 MHz					
Time Base Accuracy Limits			Measured time base error (ppm)	Pass/Fail			
Time Base Accuracy Limit: $\pm 10$ ppm							
Internal Trigger Sensitivity				Channel 1	Channel 2	Channel 3*	Channel 4*
	Generator Setting	Test Limits					
1.5 GHz models:	1.5 GHz	< 10 mV/div: greater of 1.5 div or 5 mVpp >= 10 mV/div: 1 div					
	1 GHz	< 10 mV/div: greater of 1 div or 5 mVpp >= 10 mV/div: 0.6 div					
1 GHz models:	1 GHz	< 10 mV/div: greater of 1 div or 5 mVpp					
500 MHz models:	500 MHz	>= 10 mV/div: 0.6 div					
350 MHz models:	350 MHz						
200 MHz models:	200 MHz						

External Trigger Sensitivity			
Input Range: 8 V	Generator Setting	Test Limits	Ext Trig In
	200 MHz	350 mV	_____
	100 MHz	200 mV	_____
Input Range: 1.6 V	Generator Setting	Test Limits	Ext Trig In
	200 MHz	70 mV	_____
	100 MHz	40 mV	_____
* Where applicable			





## 3 Calibrating and Adjusting

This chapter explains how to adjust the oscilloscope for optimum operating performance. You should perform self-calibration according to the following recommendations:

- Every two years or after 4000 hours of operation
- If the ambient temperature is  $>10^{\circ}\text{C}$  from the calibration temperature
- If you want to maximize the measurement accuracy

The amount of use, environmental conditions, and experience with other instruments help determine if you need shorter adjustment intervals.

### Let the Equipment Warm Up Before Adjusting

Before you start the adjustments, let the oscilloscope and test equipment warm up for at least 30 minutes.

### Read All Cautions and Warnings

Read the following cautions and warning before making adjustments or performing self-calibration.

#### **WARNING**

##### **HAZARDOUS VOLTAGES!**

Read the safety notice at the front of this book before proceeding. Maintenance is performed with power supplied to the oscilloscope and with the protective covers removed. Only trained service personnel who are aware of the hazards involved should perform the maintenance. Whenever possible, perform the procedures with the power cord removed from the oscilloscope.



#### CAUTION

##### **REMOVE POWER TO AVOID DAMAGE!**

Do not disconnect any cables or remove any assemblies with power applied to the oscilloscope. Otherwise, damage to the oscilloscope can occur.

---

#### CAUTION

##### **USE EXTERNAL FAN TO REDUCE TEMPERATURE!**

When you must operate the oscilloscope with its cover and main shield removed, use an external fan to provide continuous air flow over the samplers (the ICs with heat sinks on them). Air flow over the samplers is reduced when the cover and main shield is removed, which leads to higher than normal operating temperatures. Have the fan blow air across the system board where the heat sinks are located. If the cover is removed but the main shield remains installed and the bottom holes are not blocked, the instrument will cool properly.

---

#### CAUTION

##### **AVOID DAMAGE TO ELECTRONIC COMPONENTS!**

Electrostatic discharge (ESD) can damage electronic components. When you use any of the procedures in this chapter, use proper ESD precautions. As a minimum, place the oscilloscope on a properly grounded ESD mat and wear a properly grounded ESD strap.

---

## User Calibration

Perform user-calibration:

- Every two years or after 4000 hours of operation.
- If the ambient temperature is  $>10^{\circ}\text{C}$  from the calibration temperature.
- If you want to maximize the measurement accuracy.

The amount of use, environmental conditions, and experience with other instruments help determine if you need shorter User Cal intervals.

User Cal performs an internal self-alignment routine to optimize the signal path in the oscilloscope. The routine uses internally generated signals to optimize circuits that affect channel sensitivity, offset, and trigger parameters. Disconnect all inputs and allow the oscilloscope to warm up before performing this procedure.

Performing User Cal will invalidate your Certificate of Calibration. If NIST (National Institute of Standards and Technology) traceability is required perform the procedures in [Chapter 2](#) in this book using traceable sources.


### To perform User Cal

- 1 Disconnect all inputs from the front and rear panels, including the digital channels cable on an MSO, and allow the oscilloscope to warm up before performing this procedure.
- 2 Press the rear-panel CAL button to disable calibration protection..
- 3 Connect short (12 inch maximum) equal length cables to each analog channel's BNC connector on the front of the oscilloscope. You will need two equal-length cables for a

2-channel oscilloscope or four equal-length cables for a 4-channel oscilloscope.

Use 50 $\Omega$  RG58AU or equivalent BNC cables when performing User Cal.

- a** For a 2-channel oscilloscope, connect a BNC tee to the equal length cables. Then connect a BNC(f)-to-BNC(f) (also called a barrel connector) to the tee as shown below.

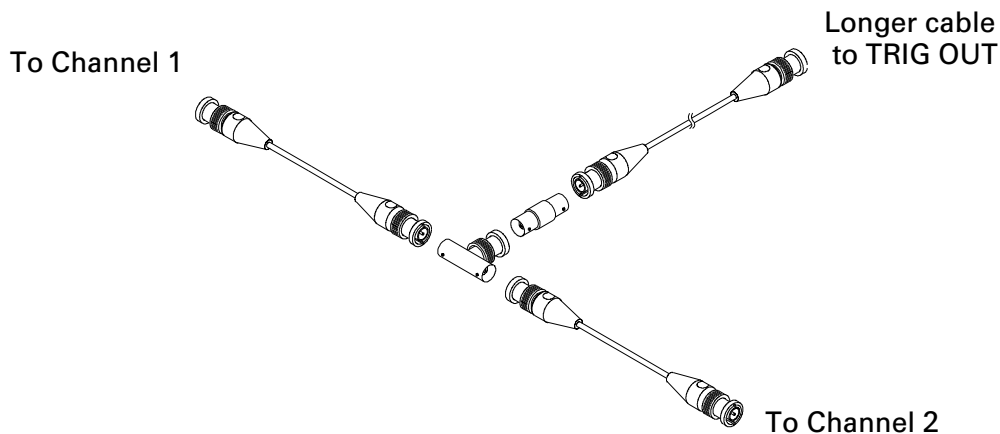


**Figure 10** User Calibration cable for 2-channel oscilloscope

- b** For a 4-channel oscilloscope, connect BNC tees to the equal-length cables as shown below. Then connect a



BNC(f)-to-BNC(f) (barrel connector) to the tee as shown below.

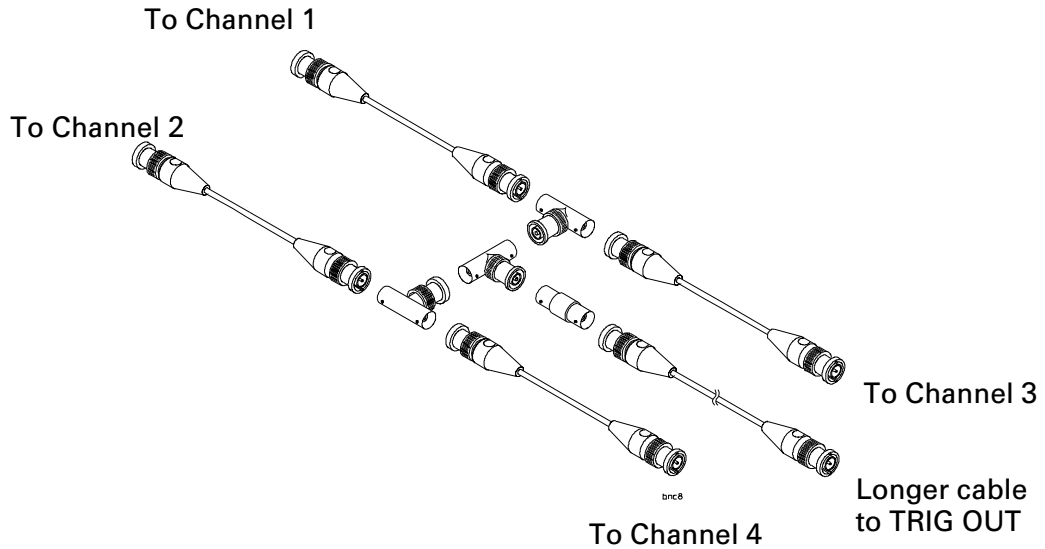

**Figure 11** User Calibration cable for 4-channel oscilloscope

- 4** Connect a BNC cable (40 inches maximum) from the TRIG OUT connector on the rear panel to the BNC barrel connector.
- 5** Press the **[Utility]** key; then, press the **Service** softkey.
- 6** Begin the Self Cal by pressing the **Start User Cal** softkey.

## **User Cal Status**

Pressing the **User Cal Status** softkey displays the following summary results of the previous User Cal, and the status of probe calibrations for probes that can be calibrated. Note that AutoProbes do not need to be calibrated, but InfiniiMax probes can be calibrated.

Results:

User Cal date:

Change in temperature since last User Cal:

Failure:

Comments:

Probe Cal Status:



## 4 Troubleshooting

Solving General Problems with the Oscilloscope 72

Verifying Basic Operation 75

Troubleshooting Internal Assemblies 82

### **Read All Cautions and Warnings**

Before you begin any troubleshooting, read all Warning and Cautions in the “Troubleshooting” section.

This chapter begins with “Solving General Problems with the Oscilloscope”. It tells you what to do in these cases:

- If there is no display.
- If there is no trace display.
- If the trace display is unusual or unexpected.
- If you cannot see a channel.

Next, this chapter describes procedures for “Verifying Basic Operation” of the oscilloscope:

- To power-on the oscilloscope.
- To perform hardware self test.
- To perform front panel self test.
- To verify default setup.
- To perform an Auto Scale on the Probe Comp signal.
- To compensate passive probes.



Finally, this chapter describes procedures for “[Troubleshooting Internal Assemblies](#)” when performing assembly-level repair:

- To prepare for internal assembly troubleshooting.
- To check the system board power supply test points.
- To check the line filter board AC output.
- To check the power switch.
- To check the power supply DC output.
- To check the display supplies.
- To check the keyboard supplies.
- To check the fan.

## Solving General Problems with the Oscilloscope

This section describes how to solve general problems that you may encounter while using the Agilent 4000 X-Series oscilloscopes.

After troubleshooting the oscilloscope, if you need to replace parts, refer to [Chapter 6](#), “Replaceable Parts,” starting on page 135.

### If there is no display

- ✓ Check that the power cord is firmly seated in the oscilloscope power receptacle.
- ✓ Check that the power source is live.
- ✓ Check that the front-panel power switch is on.
- ✓ If there is still no display, go to the troubleshooting procedures in this chapter.

## If there is no trace display

- ✓ Check that the Intensity (on the front panel) is adjusted correctly.
- ✓ Recall the default setup by pressing **[Default Setup]**. This will ensure that the trigger mode is Auto.
- ✓ Check that the probe clips are securely connected to points in the circuit under test, and that the ground is connected.
- ✓ Check that the circuit under test is powered on.
- ✓ Press the **[AutoScale]** key.
- ✓ Obtain service from Agilent Technologies, if necessary.

## If the trace display is unusual or unexpected

- ✓ Check that the Horizontal time/division setting is correct for the expected frequency range of the input signals.
- ✓ The sampling speed of the oscilloscope depends on the time/division setting. It may be that when time/division is set to slower speeds, the oscilloscope is sampling too slowly to capture all of the transitions on the waveform. Use peak detect mode.
- ✓ Check that all oscilloscope probes are connected to the correct signals in the circuit under test.
- ✓ Ensure that the probe's ground lead is securely connected to a ground point in the circuit under test. For high-speed measurements, each probe's individual ground lead should also be connected to a ground point closest to the signal point in the circuit under test.
- ✓ Check that the trigger setup is correct.
- ✓ A correct trigger setup is the most important factor in helping you capture the data you desire. See the *User's Guide* for information about triggering.

- ✓ Check that persistence in the Display menu is turned off, then press the **Clear Display** softkey.
- ✓ Press the **[Auto Scale]** key.

## If you cannot see a channel

- ✓ Recall the default setup by pressing **[Default Setup]**. This will ensure that the trigger mode is Auto.
- ✓ Check that the oscilloscope probe's BNC connector is securely attached to the oscilloscope's input connector.
- ✓ Check that the probe clips are securely connected to points in the circuit under test.
- ✓ Check that the circuit under test is powered on.

You may have pressed the **[Auto Scale]** key before an input signal was available.

Performing the checks listed here ensures that the signals from the circuit under test will be seen by the oscilloscope. Perform the remaining checks in this topic to make sure the oscilloscope channels are on, and to obtain an automatic setup.

- ✓ Check that the desired oscilloscope channels are turned on.
  - a** Press the analog channel key until it is illuminated.
  - b** On models with the MSO option, press the digital channels **[Digital]** key until it is illuminated.
- ✓ Press the **[Auto Scale]** key to automatically set up all channels.

## Verifying Basic Operation

### To power-on the oscilloscope

- 1 Connect the power cord to the rear of the oscilloscope, then to a suitable ac voltage source.

The oscilloscope power supply automatically adjusts for input line voltages in the range of 100 to 240 VAC. Ensure that you have the correct line cord (see [page 135](#)). The power cord provided is matched to the country of origin.

#### **WARNING**

##### **AVOID INJURY.**

Always operate the oscilloscope with an approved three conductor power cable. Do not negate the protective action of the three conductor power cable.

---

- Press the power switch.
    - When the oscilloscope is turned on, the front panel LEDs will briefly light up in groups from bottom to top.
    - Next the Agilent logo appears on the display.
    - Next a message will appear with tips on getting started using the oscilloscope. At this time you can press any key to remove the message and view the display. Or you can wait and the message will automatically disappear.
    - It will take a total of about 20-30 seconds for the oscilloscope to go through its basic self test and power-up routine.
- 2 Proceed to [“To perform hardware self test”](#) on page 76.

## To perform hardware self test

Pressing **[Utility] > Service > Hardware Self Test** performs a series of internal procedures to verify that the oscilloscope is operating properly.

It is recommended you run Hardware Self Test:

- After experiencing abnormal operation.
- For additional information to better describe an oscilloscope failure.
- To verify proper operation after the oscilloscope has been repaired.

Successfully passing Hardware Self Test does not guarantee 100% of the oscilloscope's functionality. Hardware Self Test is designed to provide an 80% confidence level that the oscilloscope is operating properly.

## To perform front panel self test

Pressing **[Utility] > Service > Front Panel Self Test** lets you test the front panel keys and knobs as well as the oscilloscope display.

Follow the on-screen instructions.

Failures in the front panel self test indicate problems with the keyboard, keypad, or display.



## To verify default setup

The oscilloscope is designed to turn on with the setup from the last turn on or previous setup.


To recall the default setup:

- 1 Press the **[Default Setup]** key.

This returns the oscilloscope to its default settings and places the oscilloscope in a known operating condition. The major default settings are:

- **Horizontal:**
  - main mode.
  - 100 us/div scale.
  - 0 s delay.
  - center time reference.
- **Vertical:**
  - Channel 1 on.
  - 5 V/div scale.
  - dc coupling.
  - 0 V position.
  - probe factor to 1.0 if an AutoProbe probe is not connected to the channel.
- **Trigger:**
  - Edge trigger.
  - Auto sweep mode.
  - 0 V level.
  - channel 1 source.
  - dc coupling.
  - rising edge slope.
  - 40 ns holdoff time.
- **Display:**
  - 20% grid intensity.
  - persistence off.

- **Other:**
  - Acquire mode normal.
  - Run/Stop to Run.
  - cursor measurements off.

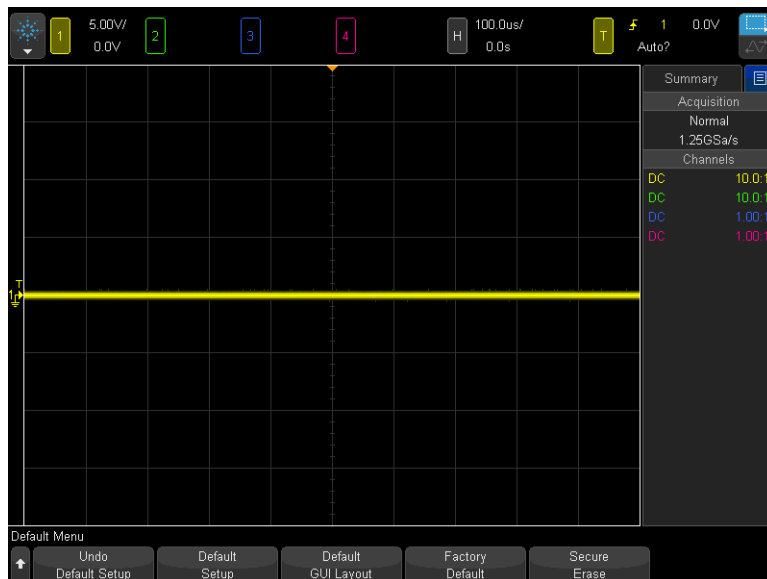


**Figure 12** Default setup screen

- 2 If your screen looks substantially different, replace the system board.

### To perform an Auto Scale on the Probe Comp signal

- 1 Press the **[Default Setup]** key. The oscilloscope is now configured to its default settings.
- 2 Connect an oscilloscope probe from channel 1 to the **Probe Comp** signal terminal on the front panel.
- 3 Connect the probe's ground lead to the ground terminal that is next to the **Demo 2** (Probe Comp) terminal.

- 4 Press **[AutoScale]**.
- 5 You should see a waveform on the oscilloscope's display similar to this:

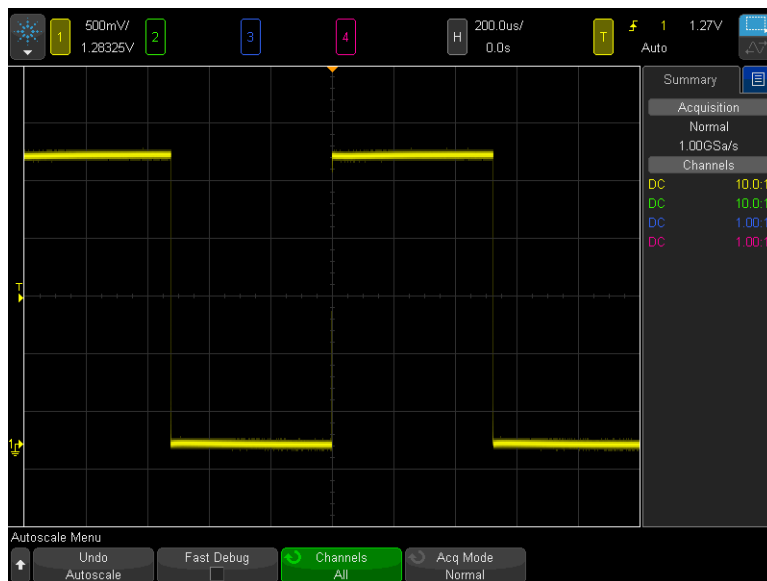

If you see the waveform, but the square wave is not shaped correctly as shown above, perform the procedure [“To compensate passive probes”](#) on page 80.

If you do not see the waveform, ensure your power source is adequate, the oscilloscope is properly powered-on, and the probe is connected securely to the front-panel analog channel input BNC and to the Demo 2 (Probe Comp) terminal.

- 6 If you still do not see the waveform, use the troubleshooting flowchart in this chapter to isolate the problem.

## To compensate passive probes

You should compensate your passive probes to match their characteristics to the oscilloscope's channels. A poorly compensated probe can introduce measurement errors.

- 1 Perform the procedure [“To perform an Auto Scale on the Probe Comp signal”](#) on page 78
- 2 Press the channel key to which the probe is connected ([1], [2], etc.).
- 3 In the Channel Menu, press **Probe**.
- 4 In the Channel Probe Menu, press **Probe Check**; then, follow the instructions on-screen.

If necessary, use a nonmetallic tool (supplied with the probe) to adjust the trimmer capacitor on the probe for the flattest pulse possible.

On the N2862/63/90 probes, the trimmer capacitor is the yellow adjustment on the probe tip. On other probes, the trimmer capacitor is located on the probe BNC connector.

**Perfectly compensated**



**Over compensated**



**Under compensated**



comp.cdr

**Figure 13** Example pulses

- 5 Connect probes to all other oscilloscope channels (channel 2 of a 2-channel oscilloscope, or channels 2, 3, and 4 of a 4-channel oscilloscope).

**6** Repeat the procedure for each channel.

The process of compensating the probes serves as a basic test to verify that the oscilloscope is functional.

# Troubleshooting Internal Assemblies

The service policy for all bandwidth model oscilloscopes is assembly level replacement. You can use the procedures described in this section to help identify assemblies that need replacement.

Generally, you want to make sure cables to the assembly are good and properly seated. Then, you check that the assembly is supplied with the proper power. If cables are good and the power is good, but the assembly still does not function properly, it must be replaced.

If you need parts or assistance from Agilent Technologies to repair your instrument, go to [www.agilent.com](http://www.agilent.com) and locate the service facility for your area.

## Equipment Required for Troubleshooting Internal Assemblies

The equipment listed in this table is required to troubleshoot the oscilloscope.

**Table 18** Equipment Required to Troubleshoot the Oscilloscope

Equipment	Critical Specifications	Recommended Model/Part
Digital multimeter	Accuracy $\pm 0.05\%$ 1 mV resolution	Agilent 34401A
Oscilloscope	Capable of measuring $\geq 500$ MHz signal. 1 M $\Omega$ input impedance.	Agilent DS06102A, MS06102A, DS07104A/B, or MS07104A/B

## To prepare for internal assembly troubleshooting

**WARNING****HAZARDOUS VOLTAGES EXIST — REMOVE POWER FIRST !**

The procedures described in this section are performed with power supplied to the oscilloscope and with the protective covers removed. Only trained service personnel who are aware of the hazards involved should perform the procedures. Whenever possible, perform the procedures with the power cord removed from the oscilloscope. Read the safety notice at the back of this book before proceeding.

---

**WARNING****HAZARDOUS VOLTAGES EXIST — HIGH VOLTAGE IS PRESENT ON POWER SUPPLY HEAT SINKS !**

The power supply heat sinks of the 4000 X-Series oscilloscopes are at a high potential. This presents an electric shock hazard. Protect yourself from electric shock by keeping this area covered or by not coming in contact with the heat sinks when the power cord is attached to the oscilloscope!

---

**CAUTION****REMOVE POWER TO AVOID DAMAGE !**

Do not disconnect any cables or remove any assemblies while power is applied to the oscilloscope, or damage to the oscilloscope can occur.

---


**CAUTION****AVOID ESD DAMAGE TO COMPONENTS !**

ELECTROSTATIC DISCHARGE (ESD) can damage electronic components. Use proper ESD precautions when doing any of the procedures in this chapter. As a minimum, place the oscilloscope on a properly grounded ESD mat and wear a properly grounded ESD strap.

---

- 1 Disconnect any external cables from the front panel.
- 2 Disconnect the power cord.
- 3 Remove the cabinet following the instructions on [page 104](#).
- 4 Remove the power supply shield and air duct following the instructions on [page 126](#).
- 5 Separate the front and rear decks following the instructions on [page 105](#), but leave all cables except the

fan power cable connected as shown in the following figure.

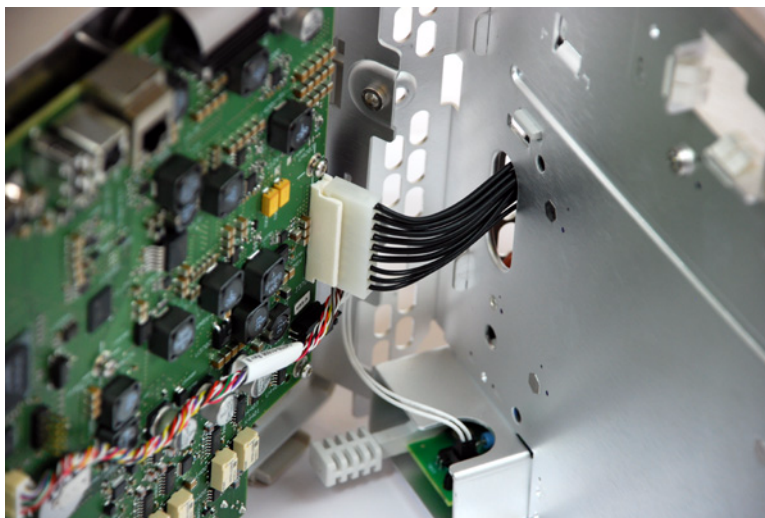

**Figure 14** Setup for troubleshooting internal assemblies

**CAUTION**

**MAKE SURE EARTH GROUND IS MAINTAINED FOR THE FRONT DECK !**

For example, connect a cable with alligator clips between the rear deck chassis and the front deck chassis, or connect a BNC cable from one of the BNCs on the system board to a known grounded BNC on your workbench.

Other advice:

- Place the front deck in a cover or rest it on a cloth to prevent scuffing the front panel knobs.
- Because of the short length of the interboard supply cable, elevate the rear deck to prevent strain on the cable.

- 6 Make sure the keyboard cable, display cable, display backlight power cable, and all other cables except the fan power cable are properly connected.



**CAUTION****USE AN EXTERNAL FAN TO AVOID OVERHEATING COMPONENTS !**


When you remove the oscilloscope cover and main shield, use an external fan to provide continuous air flow over the heat sinks. Air flow over the heat sinks is reduced when the cover and main shield are removed, which leads to higher than normal operating temperatures. Have the fan blow air across the system board where the heat sinks are located. Otherwise, damage to the components can occur.

If the cabinet of a 4000 X-Series oscilloscope is removed but the main power supply shield remains installed and the bottom holes are not blocked, the instrument will cool properly.

---


# Flowchart for Internal Assembly Power Troubleshooting

The following flowchart is a simplified overview of troubleshooting power to the oscilloscope's internal assemblies.

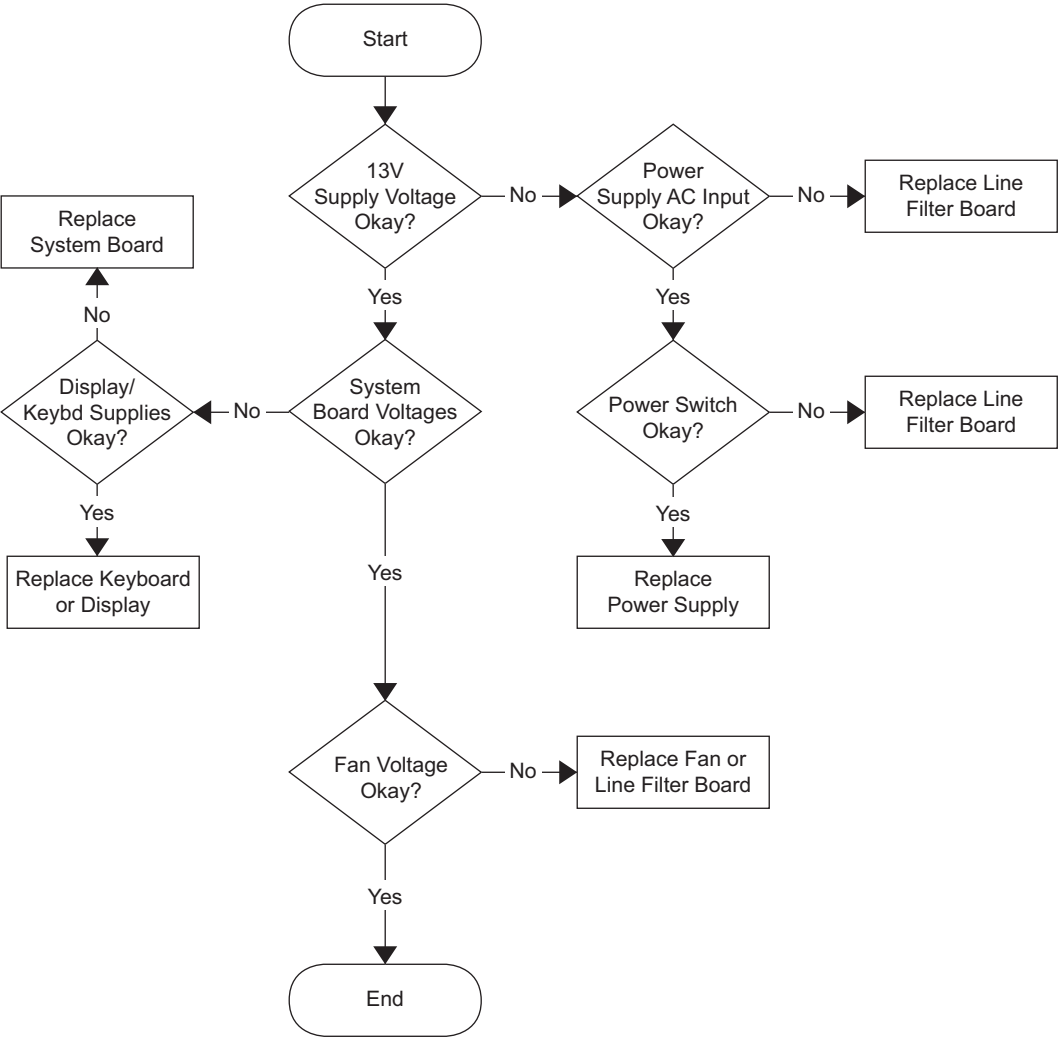


## System Board Drawings

Use these drawings to locate test points on the oscilloscope system board.

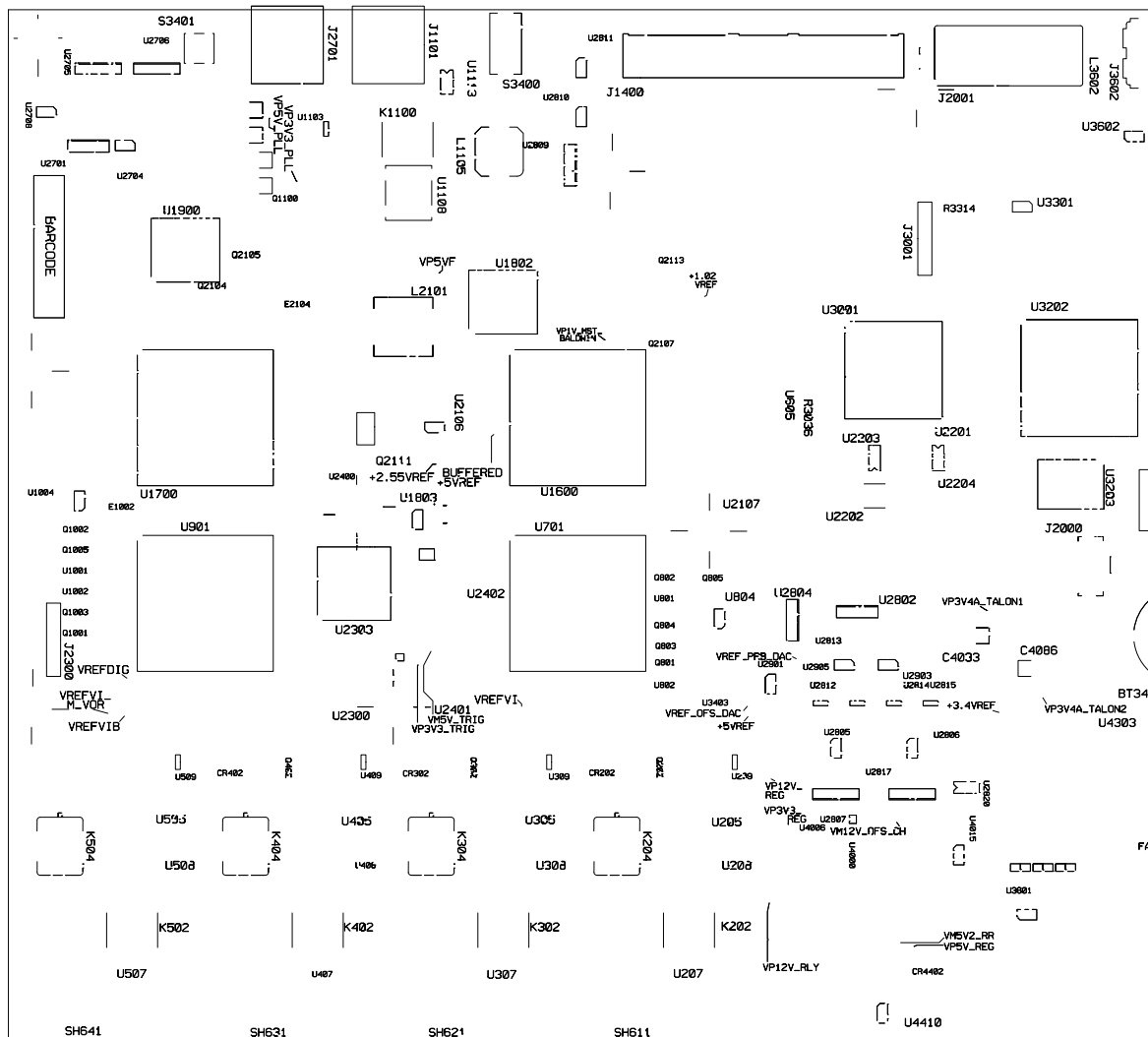

**Figure 15** System Board Test Points/Connectors - Top Side, Left

## 4

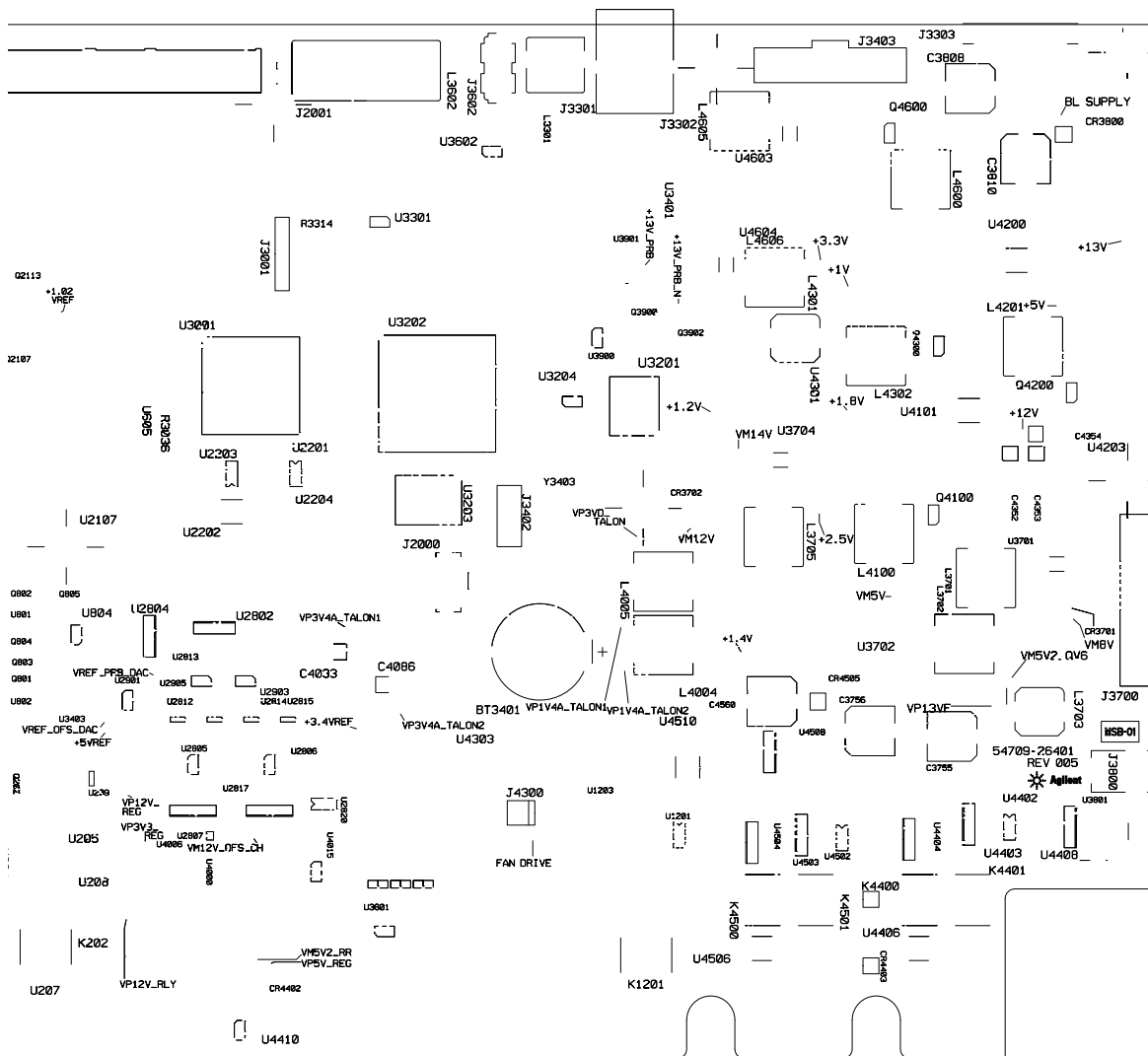

**Figure 16** System Board Test Points/Connectors - Top Side, Right

## To check the system board power supply test points

This procedure checks the power supply test points on the system board (see See [Figure 16](#) on page 88). Values outside the expected range help identify bad assemblies.

- 1 Follow the instructions in [“To prepare for internal assembly troubleshooting”](#) on page 83.
- 2 Connect the negative lead of the multimeter to a ground point on the oscilloscope.
- 3 Connect the power cord, and turn on the oscilloscope.
- 4 First check the bulk power supply voltage:

Test Point (near J3700)	Expected Value	Assemblies Supplied	If Good	If Bad
VP13V (+13V)	+13 V $\pm$ 3%	All	Indicates power supply is good.	Go to <a href="#">“To check the line filter board AC output”</a> on page 91.


5 Next, check the supplies coming from the system board:

Test Point (near J3700)	Expected Value	Assemblies Supplied	If Good	If Bad
VP5V (+5V)	+5 V $\pm 3\%$	Keyboard, System Board	Indicates supply coming from system board is good.	Replace the system board.
VP12V (+12V)	+12 V $\pm 3\%$	System Board		Replace the system board.
VP1V8 (+1.8V)	+1.8 V $\pm 3\%$	System Board		Replace the system board.
VP1V4 (+1.4V)	+1.4 V $\pm 3\%$	System Board		Replace the system board.
VP1V (+1V)	+1 V $\pm 3\%$	System Board		Replace the system board.
VP1V2 (+1.2V)	+1.2 V $\pm 3\%$	System Board		Replace the system board.
VP2V5 (+2.5V)	+2.5 V $\pm 3\%$	System Board		Replace the system board.
VP3V3 (+3.3V)	+3.3 V $\pm 3\%$	LCD, System Board		Go to <a href="#">“To check the display supplies”</a> on page 95.
VM12V	-11.5 V to -12.5 V	System Board		Replace the system board.
VM14V	-13.5 V to -14.5 V	System Board		Replace the system board.
VM8V	-7.814 V to -8.333 V	System Board		Replace the system board.
VM5V2	-5.023 V to -5.329 V	System Board		Replace the system board.
VP13VF	+13 V $\pm 3\%$	System Board		Replace the system board.
VP25V (BL SUPPLY)	+24.84 V to +25.2 V	Display Backlight Power		Go to <a href="#">“To check the display supplies”</a> on page 95.
VP13V_PRB (+13V_PRB)	+13 V $\pm 3\%$	AutoProbe Interface		Go to <a href="#">“To check the keyboard supplies”</a> on page 96.

## To check the line filter board AC output

When the 13 V bulk power is not being properly supplied, this procedure tests the AC input to the power supply to determine whether there is a problem with the line filter board.

- 1 Follow the instructions in [“To prepare for internal assembly troubleshooting”](#) on page 83.
- 2 Connect the power cord, and turn on the oscilloscope.
- 3 Verify that AC power is present at J1 on the power supply using a DVM and probes as shown in the following picture (remember this is an AC voltage measurement!).

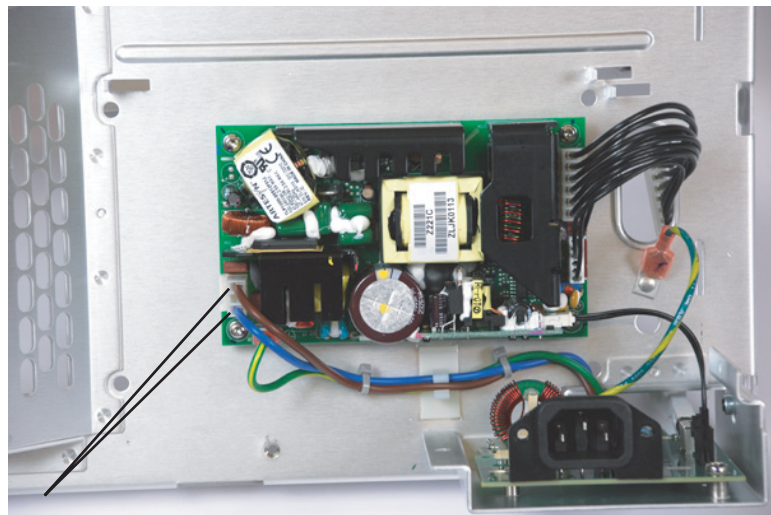

Verify AC potential between these two points

**Figure 17** Verify line filter board AC output

- If you have AC power equal to what is being applied to J101 (power cord socket) on the line filter assembly, the AC mains portion of the line filter assembly is probably okay.


- If there is no AC power at J1 of the power supply, there is something wrong with the AC mains section of the line filter assembly, and you need to *replace the line filter assembly*.
- If there is AC power at J1 of the power supply, but the instrument still will not power ON when the power switch (S101) is in the ON position, go to [“To check the power switch”](#) on page 92.

### To check the power switch

This procedure verifies the operation of the power switch (S101) on the line filter board assembly.

- 1 Follow the instructions in [“To prepare for internal assembly troubleshooting”](#) on page 83.
- 2 With the power cord disconnected, verify the operation of the power switch (S101) using a DMM in the resistance measurement mode and a set of probes.
  - a Remove the “remote on/off” cable from connector J107 on the line filter board.
  - b Place the DMM probes on pins 1 and 2 of J107.



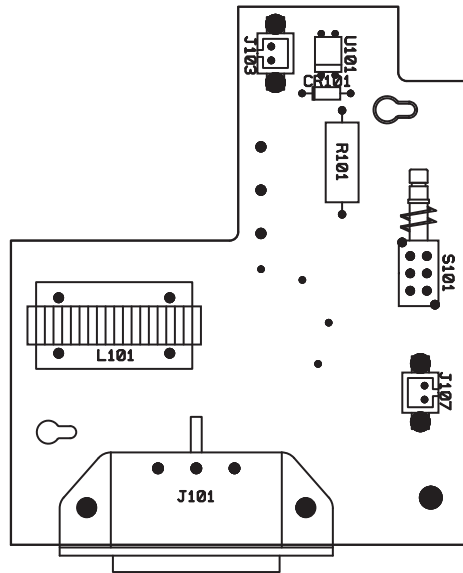

**Figure 18** Verify power switch operation

- With S101 in the OFF position (switch contacts closed) you should measure less than 1 ohm of resistance.
- With S101 in the ON position (switch contacts open) you should measure infinite resistance.

If you do not measure these two values, there is something wrong with the power switch (S101 on the line filter assembly), and you need to *replace the line filter assembly*.

If the power switch is good, and the AC input to the power supply is good, but the 13 V bulk power is not being properly supplied, it is likely that the power supply is bad. To rule out a bad DC supply cable, see [“To check the power supply DC output”](#) on page 94.

## To check the power supply DC output

This procedure checks the power supply DC output after verifying the AC input is good and checking the power switch operation.

- 1 Follow the instructions in [“To prepare for internal assembly troubleshooting”](#) on page 83.
- 2 Connect the negative lead of the multimeter to a ground point on the oscilloscope.
- 3 With the power cord disconnected:
  - a Disconnect the DC supply cable that is plugged into J3700 of the system board.
- 4 Connect the power cord, and turn on the oscilloscope.
- 5 Verify there is +13V between pins 1 (-) and 6 (+) of J3 on the power supply.

If +13V is not present between pins 1 and 6 of J3 on the power supply, *replace the power supply*.

If +13V is present between pins 1 and 6 of J3 on the power supply, either the cable is wired incorrectly or the power switch S101 is defective or has a solder short to ground between one or some of its pins.

- a Turn off the oscilloscope, and disconnect the power cord.
- b Replace the DC supply cable.
- c Connect the power cord, and turn on the oscilloscope.
- d Check again to see if the output of the power supply turns ON and OFF as you actuate the power switch S101.

If replacing the cable does not fix the problem, there is something wrong with the power switch S101, and you must *replace the line filter assembly*.

## To check the display supplies

This procedure checks the display supply voltages on the system board test points when the display backlight power cable and the display cable are disconnected.

If one or both of these voltages are bad when the cables are connected but good when the cables are disconnected, it indicates problems with the display assembly.

- 1 Follow the instructions in [“To prepare for internal assembly troubleshooting”](#) on page 83.
- 2 Connect the negative lead of the multimeter to a ground point on the oscilloscope.
- 3 With the power cord disconnected:
  - a Disconnect the display backlight power cable from the system board connector J3303.
  - b Disconnect the display cable from the system board connector J2000.
- 4 Connect the power cord, and turn on the oscilloscope.
- 5 Check the display supplies:

Test Point (near J3303)	Expected Value	Assemblies Supplied	If Good	If Bad
VP25V (BL SUPPLY)	+24.84 V to +25.2 V	Display Backlight Power	Replace the display assembly.	Replace the system board.

Test Point (near J2000)	Expected Value	Assemblies Supplied	If Good	If Bad
VP3V3 (+3.3V)	+3.3 V $\pm$ 3%	LCD, System Board	If good when the display cable is disconnected from the system board, but bad when the display cable is connected, this indicates a problem with either the display cable or the display assembly. <b>1</b> Turn off the oscilloscope, and disconnect the power cord. <b>2</b> Replace the display cable. <b>3</b> Connect the power cord, and turn on the oscilloscope. If the problem is not fixed, replace the display assembly.	Replace the system board.

## To check the keyboard supplies

This procedure checks the keyboard supply voltages on the system board test points when the keyboard cable is disconnected.

If one or both of these voltages are bad when the keyboard cable is connected but good when the cable is disconnected, it indicates problems with the keyboard assembly.

- 1** Follow the instructions in [“To prepare for internal assembly troubleshooting”](#) on page 83.
- 2** Connect the negative lead of the multimeter to a ground point on the oscilloscope.
- 3** With the power cord disconnected:
  - a** Disconnect the keyboard cable from the system board connector J3403.
- 4** Connect the power cord, and turn on the oscilloscope.

5 Check the keyboard supplies:

Test Point (near J4303)	Expected Value	Assemblies Supplied	If Good	If Bad
VP5V (+5V)	+5 V $\pm 3\%$	Keyboard, System Board	If good when the keyboard cable is disconnected from the system board, but bad when the keyboard cable is connected, this indicates a problem with either the keyboard cable or the keyboard,	Replace the system board.
VP13V_PRB (+13V_PRB)	+13 V $\pm 3\%$	AutoProbe Interface		Replace the system board.
VM15V_PRB (-15V_PRB)	-13.2 V $\pm 4\%$	AutoProbe Interface	1 Turn off the oscilloscope, and disconnect the power cord.	
VP7V5_PRB (+7.5_PRB)	+6.28 V $\pm 4\%$	AutoProbe Interface	2 Replace the keyboard cable. 3 Connect the power cord, and turn on the oscilloscope.	
VM7V5_PRB (-7.5_PRB)	-6.25 V $\pm 4\%$	AutoProbe Interface	a If the problem is not fixed, replace the keyboard.	

### To check the fan


The fan speed is controlled by a circuit on the system board.

If the fan is running, perform the hardware self-tests. Go to [“To perform hardware self test”](#) on page 76.

If the fan is not running, it may be defective. Follow these steps:

- 1 Follow the instructions in [“To prepare for internal assembly troubleshooting”](#) on page 83.
- 2 Disconnect the fan cable from the system board.
- 3 Connect the power cord, and turn on the oscilloscope.
- 4 Measure the fan voltage at the J4300 connector on the system board.


See the following figure for the location of the fan connector.

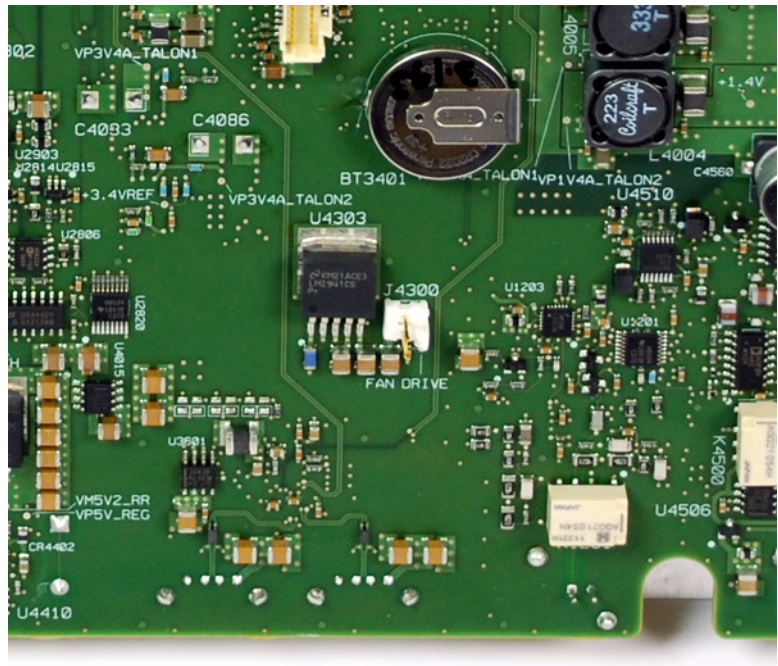

- 5 If the fan voltage is approximately +9 Vdc at room temperature, replace the fan. If the fan voltage is not approximately +9 Vdc, *replace the system board.*

The proper voltage range depending on temperature is between +7.9 Vdc to +12.8 Vdc.

**Figure 19** Location of the Fan Connector

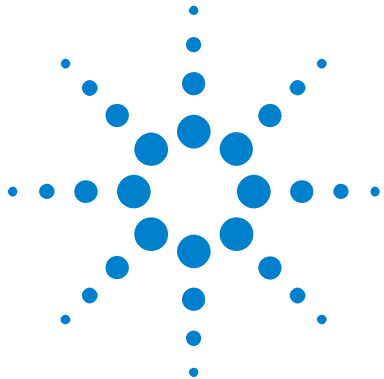

## 5 Replacing Assemblies

The service policy for 200 MHz and lower bandwidth oscilloscopes is unit replacement, so the instructions in this chapter are primarily for the 350 MHz, 500 MHz, and 1 GHz bandwidth oscilloscopes.

This chapter describes how to remove assemblies from an oscilloscope. To install a replacement assembly after you have removed an old one, follow the instructions in reverse order.

The parts shown in the following figures are representative and may look different from what you have in your oscilloscope.

### **Instructions for removable assemblies include:**

- To remove the handle [102](#)
- To remove the adjustable legs [103](#)
- To remove the cabinet [104](#)
- To remove the rear deck assembly [105](#)
- To replace the battery [108](#)
- To remove the acquisition board [109](#)
- To remove the touch controller board [112](#)
- To remove the front panel knobs [115](#)
- To remove the front bezel assembly [116](#)
- To remove the display assembly [120](#)
- To remove the keyboard and keypad [122](#)
- To remove the fan assembly [125](#)
- To remove the power supply shield [126](#)



To remove the line filter board [129](#)

To remove the power supply [132](#)

### Tools Used for Disassembly

Use these tools to remove and replace the oscilloscope assemblies:

- T6, T10, and T20 TORX drivers.
- 5/8-inch socket driver or adjustable wrench (for BNC nuts).
- Flat head screw driver.

### See how the Oscilloscope Parts Fit Together

An exploded view of the oscilloscope is included in the “Replaceable Parts” chapter. It shows the individual part numbers used in the assemblies, and shows you how the parts fit together.

### Read All Warnings and Cautions

Read the following warnings and cautions before removing and replacing any assemblies in the oscilloscope.

#### WARNING

##### HAZARDOUS VOLTAGES!

Read the safety summary at the back of this book before proceeding. Maintenance is performed with power supplied to the oscilloscope and with the protective covers removed. Only trained service personnel who are aware of the hazards involved should perform the maintenance. Whenever possible, perform the procedures with the power cord removed from the oscilloscope.

#### WARNING



##### AVOID ELECTRICAL SHOCK!

Hazardous voltages exist on the LCD assembly and power supply. To avoid electrical shock:

- 1 Disconnect the power cord from the oscilloscope.
- 2 Wait at least three minutes for the capacitors in the oscilloscope to discharge before you begin disassembly.

Read the Safety Summary at the back of this manual before you begin.



**CAUTION****REMOVE POWER TO AVOID DAMAGE!**

Remove power before you begin to remove and replace assemblies. Do not remove or replace assemblies while the oscilloscope is turned on, or damage to the components can occur.

---

**CAUTION****AVOID DAMAGE TO ELECTRONIC COMPONENTS!**

ELECTROSTATIC DISCHARGE (ESD) can damage electronic components. When doing any of the procedures in this chapter, use proper ESD precautions. As a minimum, you should place the instrument on a properly grounded ESD mat and wear a properly grounded ESD strap.

---

## To remove the handle

The strap handle must be removed prior to removing the cabinet. The removal of the Strap handle may also be necessary when mounting oscilloscope in a rack.

- 1 Using T20 TORX driver, remove the screws holding strap handle caps in place.
- 2 Lift the strap handle off the cabinet.



**Figure 20** Removing the handle

## To remove the adjustable legs

Adjustable legs must be removed prior to removing the cabinet.

- 1 Using a T20 TORX, remove shoulder screw and washer.
- 2 Pull adjustable leg assembly from cabinet. Assembly contains latching mechanism and spring.




**Figure 21** Removing adjustable legs

## To remove the cabinet

Removing the cabinet allows access to the rear deck, fan assembly, power supply cover, power supply assembly, and power switch assembly.

- 1 Using T20 TORX, remove the three screws securing cabinet to rear deck assembly.
- 2 Carefully slide cabinet back away from rear deck assembly.

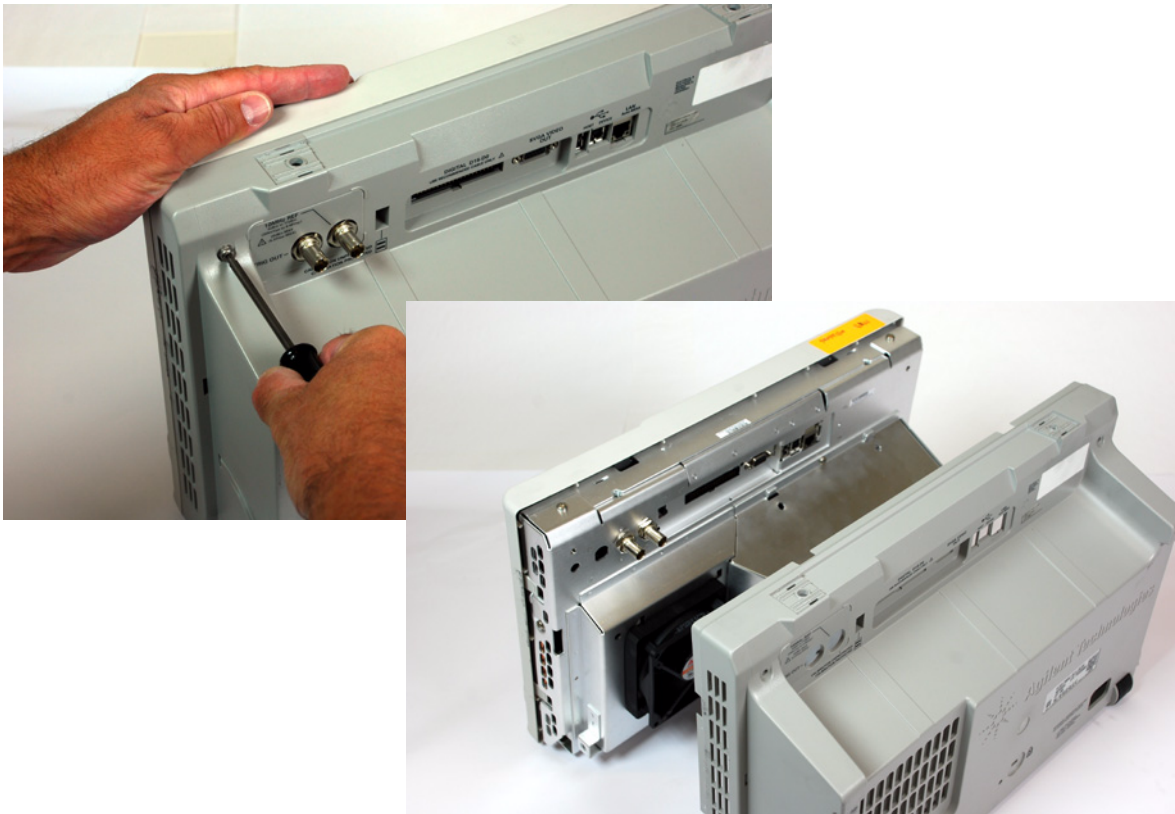

**Figure 22** Removing the cabinet

## To remove the rear deck assembly

Removing the rear deck allows access to the front deck, acquisition board, and inverter board.





- 1 Remove the BNC securing nuts and washers.



**Figure 23** Removing the BNC securing nuts and washers

## 5 Replacing Assemblies

- 2 Using a T20 TORX, locate and remove all screws securing rear deck to front deck.

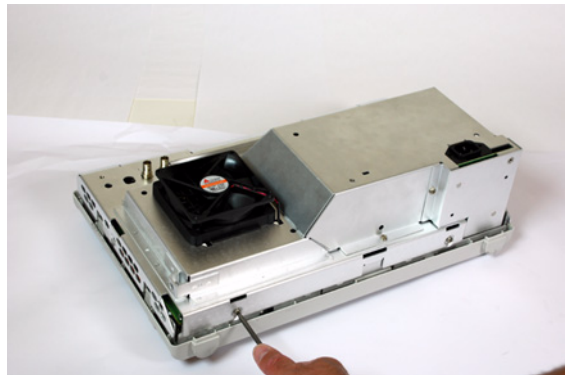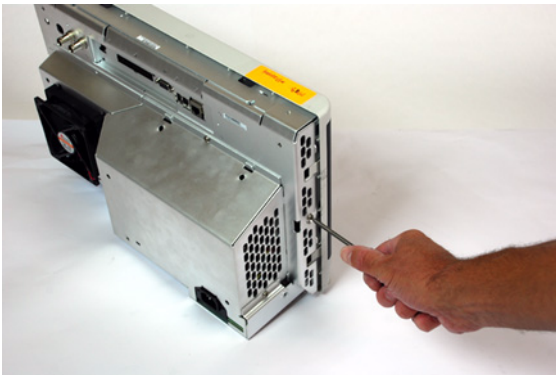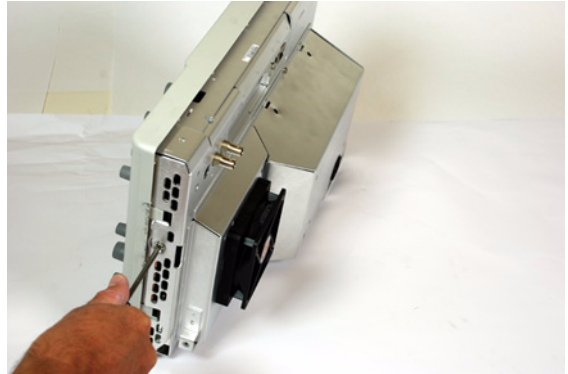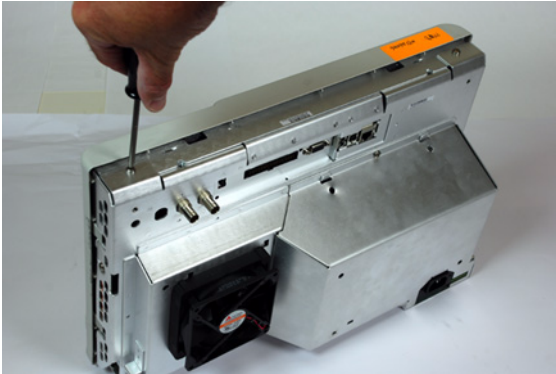

**Figure 24** Removing the rear deck assembly


- 3 Carefully separate rear deck from front deck. Take care not damage extender switch.

### **WARNING**

Sheet metal parts may have sharp edges. Handle with care to avoid injury.




- 4 Disconnect power harness and fan cable from acquisition board. Note cable locations for re-assembly

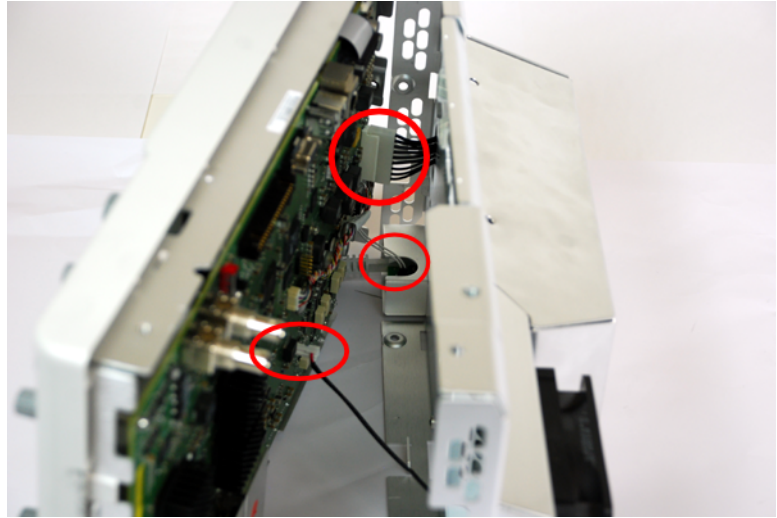

**Figure 25** Separating front and rear deck assemblies

## To replace the battery

If the battery needs to be replaced, use a CR2032/1HG or CR2032/HGN 3V manganese dioxide lithium battery.

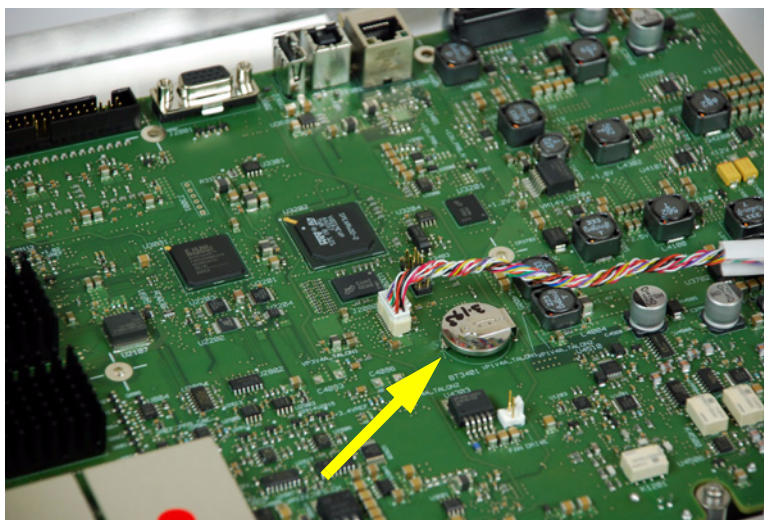

**Figure 26** Battery location



## To remove the acquisition board

The following illustrates how to remove the Acquisition and Dual inverter printed circuit boards.

- 1 Using a TORX T6 driver locate and remove the 4 screws on the front of the instrument (4 Channel version).




**Figure 27** Removing the analog channel BNC securing T6 screws

- 2 Disconnect Inverter, keyboard and display cables. Note locations for re-connection. It should be noted that cables can be removed from cable clamps at this time as well. The cable restraining pads that affix the display cable to the front deck are adhesive and great care should be taken when removing them so as not to damage the cable.


## 5 Replacing Assemblies

- 3 Using a TORX T10 driver locate and remove 15 mounting screws.

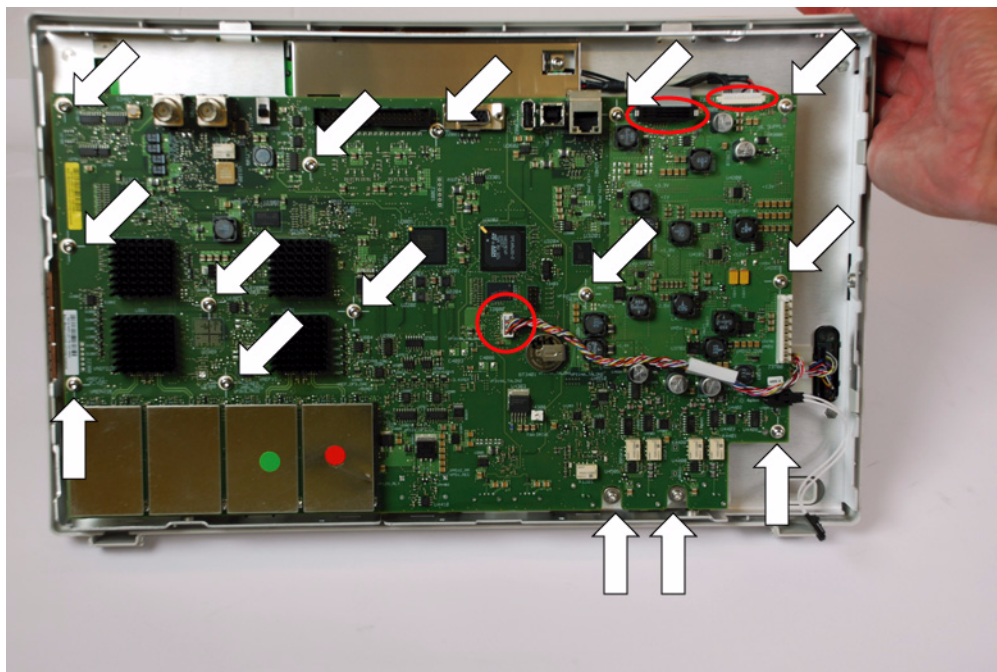

**Figure 28** Acquisition board TORX T10 mounting screw locations


- 4 Carefully lift acquisition board off front deck.

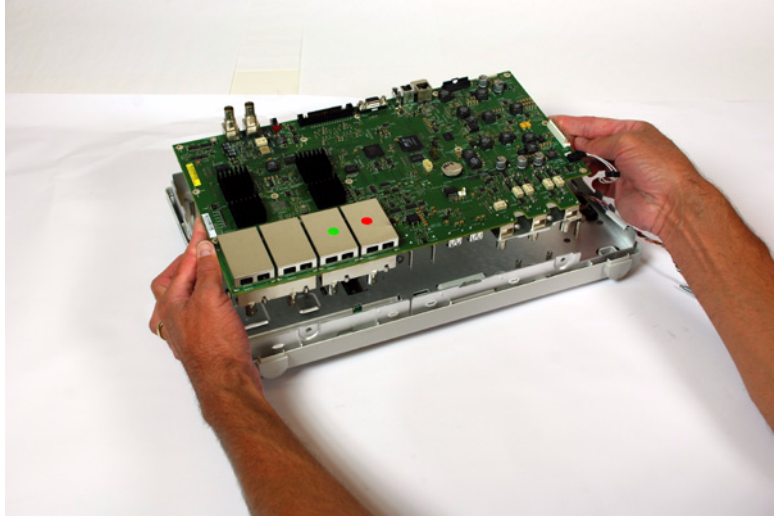


**Figure 29** Lifting acquisition board off the front deck

## To remove the touch controller board

- 1 Using a TORX T10 driver locate and remove the two screws securing the touch controller shield.

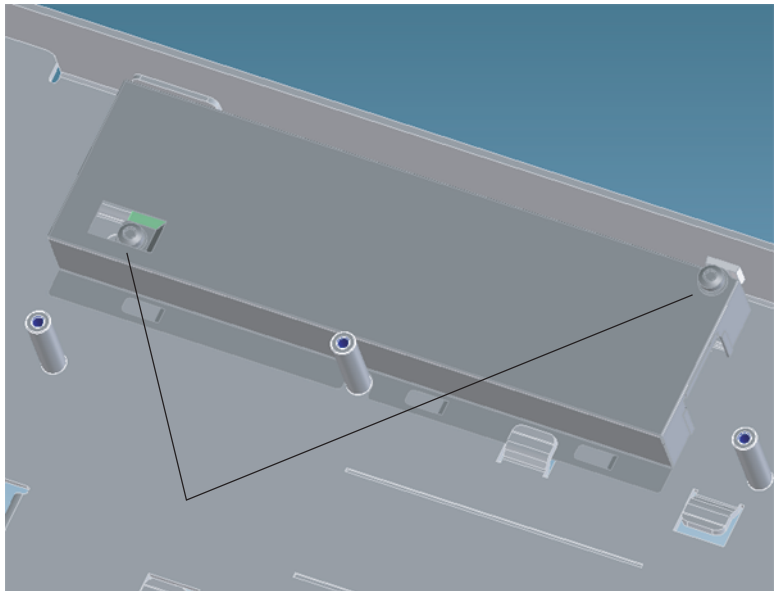


**Figure 30** Removing the touch controller shield

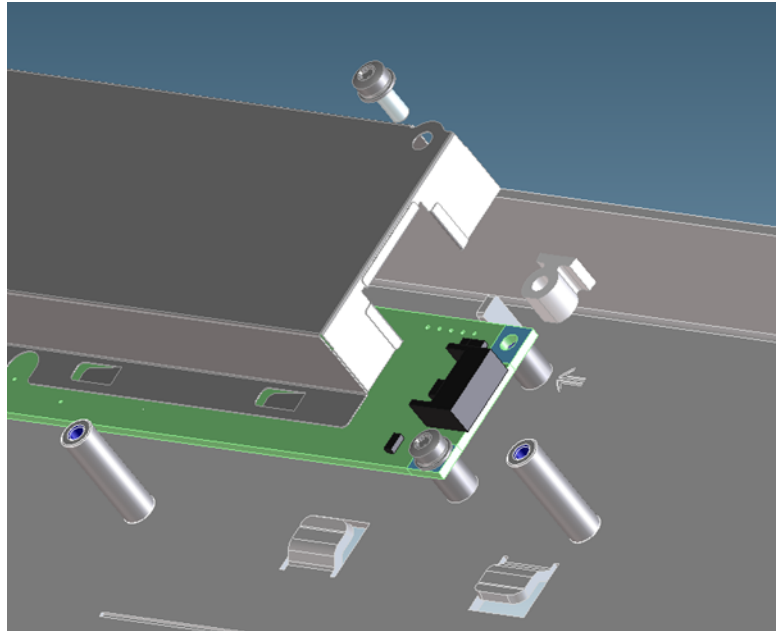

**Figure 31** Lifting the touch controller shield off

- 2 Remove the three screws securing the touch controller board to the front deck.

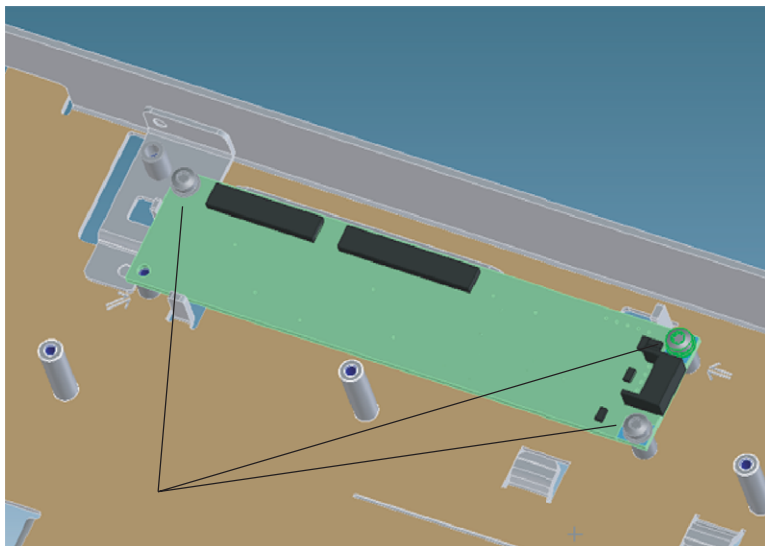

**Figure 32** Removing the touch controller board

- 3 Disconnect all cables, lift board off front deck.

## To remove the front panel knobs

The front panel knobs must be removed prior to localized front panel overlay installation or front panel disassembly.

- 1 Gently pull on the front panel knobs to remove them.



**Figure 33** Removing the front panel knobs

## To remove the front bezel assembly

- 1 Remove cable shield from front deck. Carefully squeeze the shield so that it clears the sheet metal tabs holding it place then slide it free of the front deck:

### WARNING

Thin sheet metal parts may have sharp edges. Handle with care to avoid injury.

---



**Figure 34** Removing the keyboard cable shield




- 2 Disconnect the end of the keyboard cable that extends thru front deck.



**Figure 35** Disconnecting the keyboard cable


- 3 Remove the display cable from the cable guide.

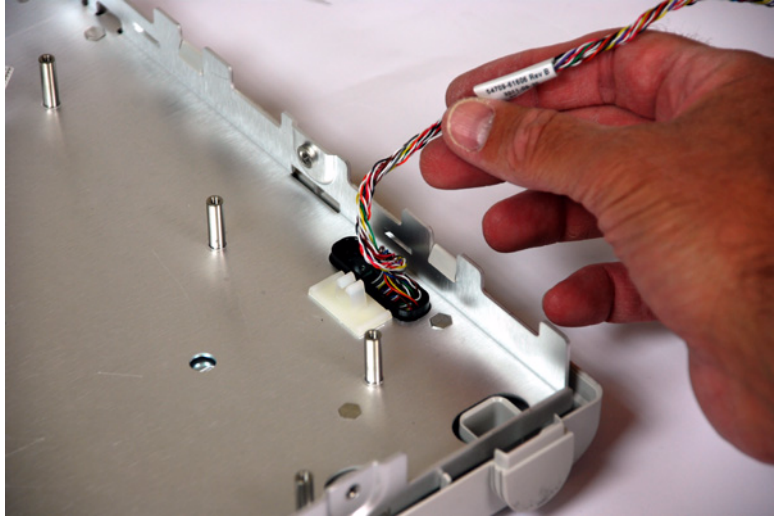

**Figure 36** Removing the display cable from the guide

**4** Remove the bezel from the front deck.


The bezel is secured to front deck by molded-in retaining clips located around the perimeter of the bezel. Gently pry these outward (either by hand or using a flat head screwdriver). Working your way around the bezel releasing the clips gently lift the bezel away from the front deck.

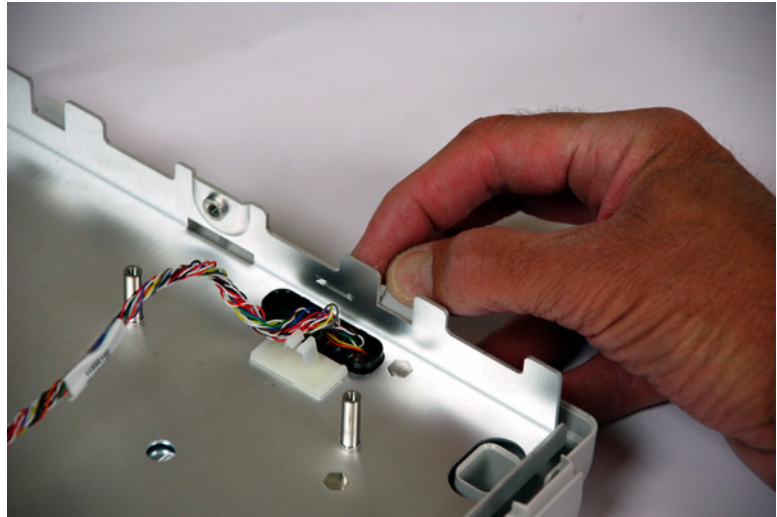

**Figure 37** Removing the bezel

## To remove the display assembly

- 1 Using a TORX T10 remove the two screws that secure the display assembly to the front deck.

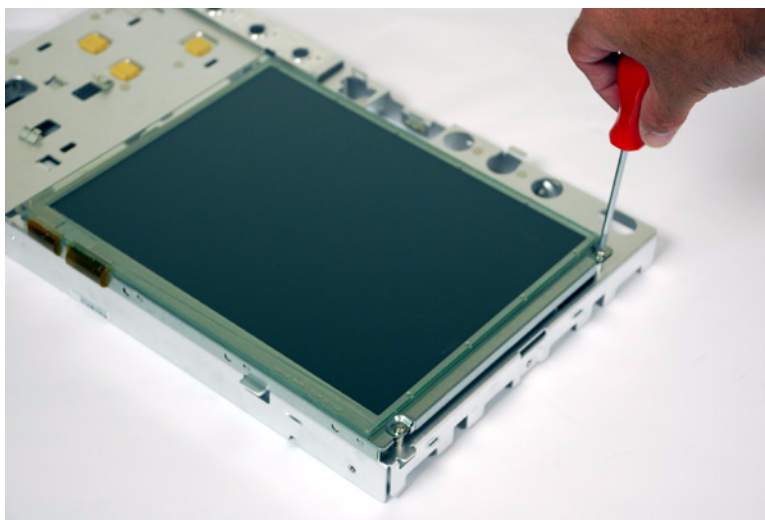

**Figure 38**    Removing the display assembly

**2** Disconnect the backlight power cable.



**Figure 39** Disconnecting the backlight power cable

Take care while threading display cables thru front deck sheet metal openings.


## To remove the keyboard and keypad

- 1 Remove the softkey board cable from the keyboard board.  
Slide the connector clamp toward the softkey board; then, slide the ribbon cable out.



**Figure 40** Disconnecting the softkey board cable

- 2 To separate the keyboard board from the bezel, carefully pull back locking tabs that secure it.

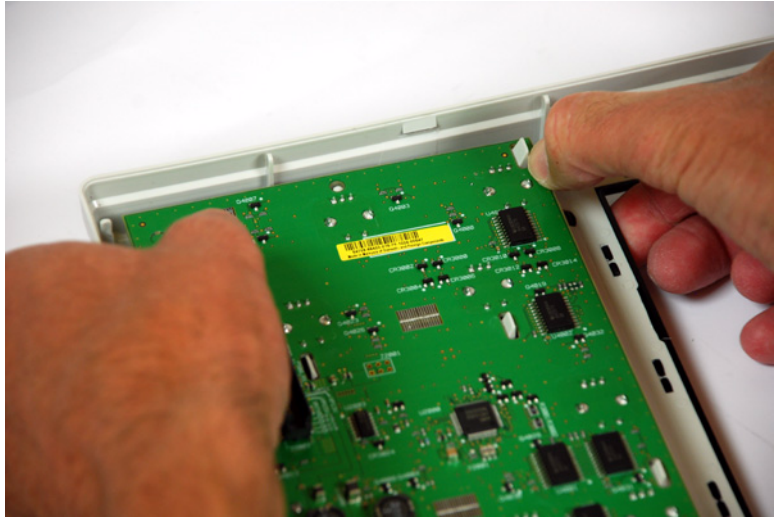


**Figure 41** Removing the keyboard




**Figure 42** Keyboard and keypad removed



- 3 To separate the softkey board from the bezel, carefully pull back locking tabs that secure it.

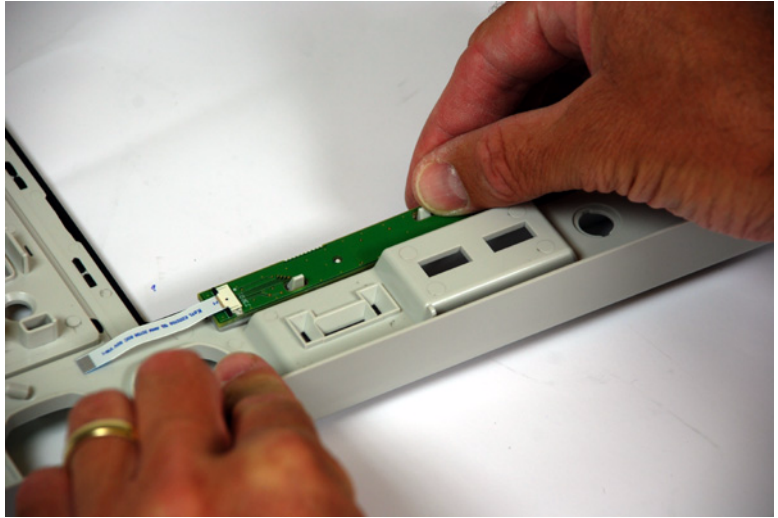

**Figure 43** Unclip the softkey board from the bezel

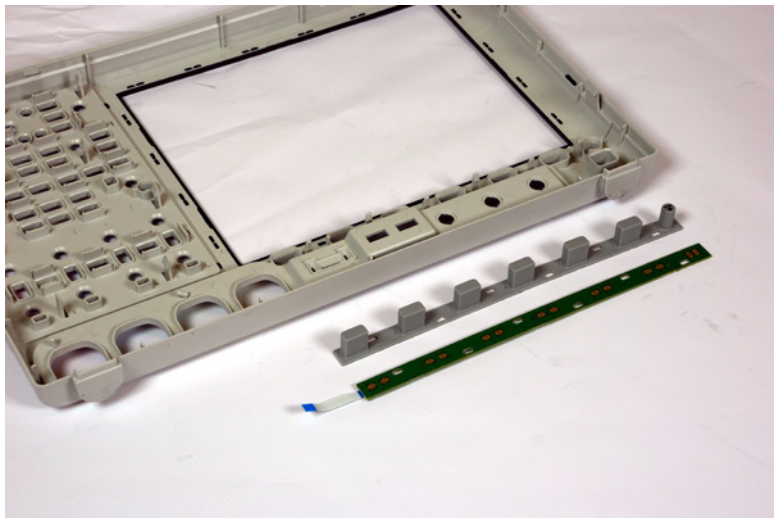


**Figure 44** Softkey board and keypad removed



## To remove the fan assembly

- 1 Ensure that fan power cable has been disconnected from acquisition board.
- 2 Carefully slide fan assembly (fan and fan mount) to the right, then lift away from rear deck.




Note, fan mount is soft and can be damaged by sharp sheet metal edges. Take care that the fan power cable is not damaged when pulling across sheet metal edges.

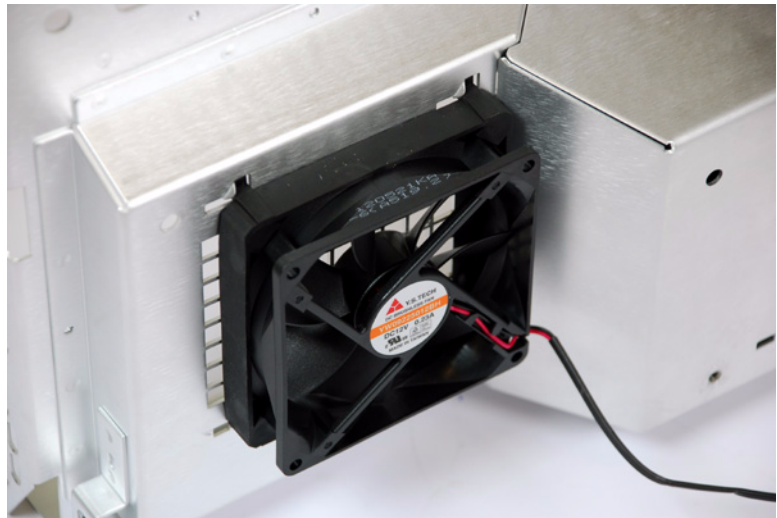

**Figure 45** Removing fan assembly

## To remove the power supply shield

- 1 To remove power supply shield, locate and remove using a TORX T20 the four screws securing the power supply shield to the rear deck.

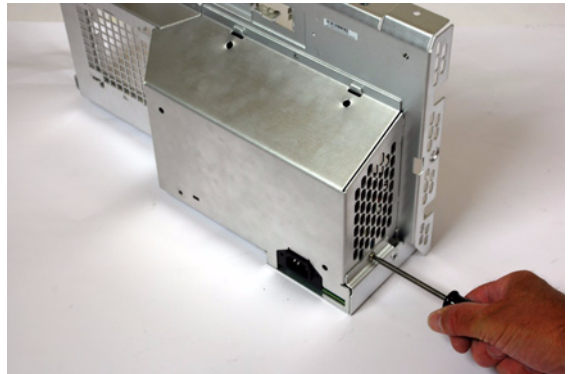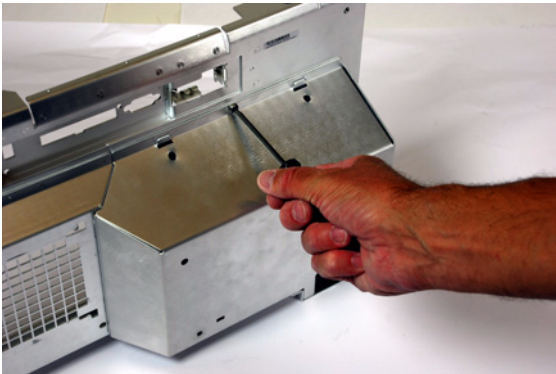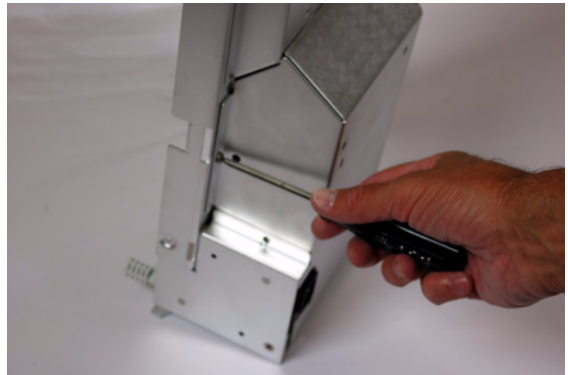


**Figure 46** Removing the power supply shield

**WARNING**

Thin sheet metal parts may have sharp edges. Handle with care to avoid injury.

---

- 2 Once screws have been removed, carefully remove the power supply cover by lifting the cover up and off retaining tabs on rear deck.

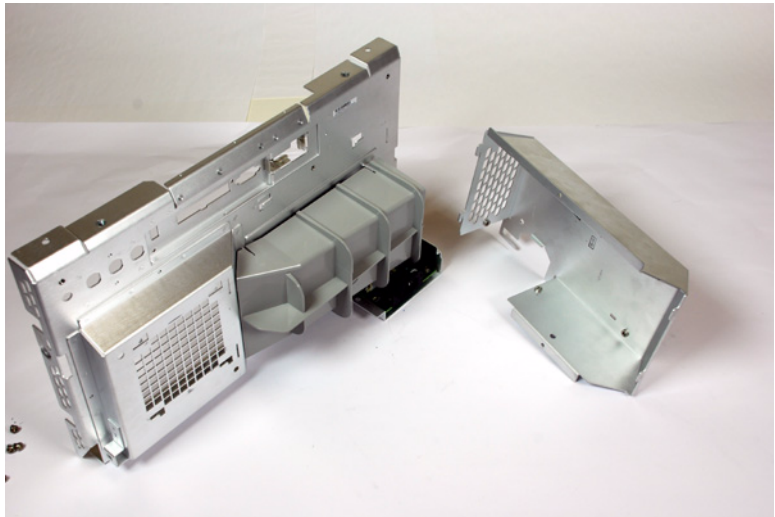

**Figure 47** Power supply shield removed


- 3 Carefully remove the air duct by lifting it up and off the rear deck.



**Figure 48** Removing the air duct

## To remove the line filter board

- 1 Disconnect the ground wire from its chassis terminal.

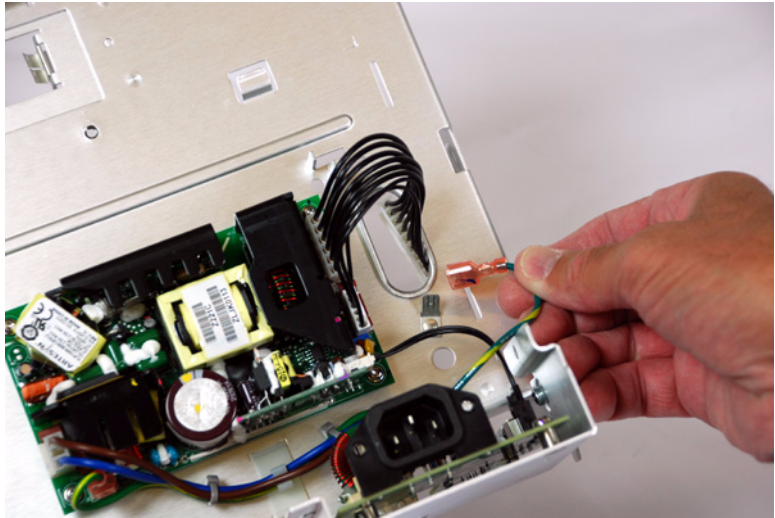

**Figure 49** Disconnecting the ground wire

- 2 Disconnect all cables from the line filter board to the power supply and remove the cables from the cable guide.

- 3 Locate and remove using a TORX T20 driver the single screw securing the assembly to the rear deck.




**Figure 50** Removing the line filter board

- 4 Slide assembly to right and lift out of rear deck.
- 5 Take care that you do not damage the switch extender during removal.

- 6 To remove the switch extender, gently pry open the extender using a flat head screwdriver.

**CAUTION**

Twisting the latch too much could cause it to break!

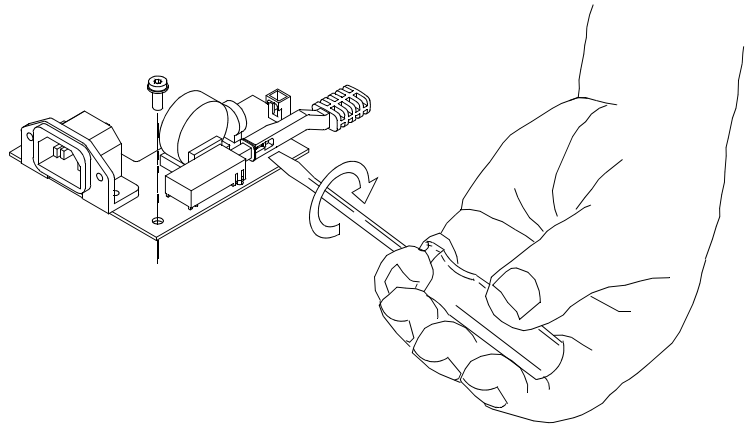

**Figure 51** Removing power switch extender

## To remove the power supply

- 1 Disconnect all cables from power supply board.
- 2 Locate and remove using a TORX T10 driver the four screws securing the power supply assembly to the rear deck.




**Figure 52** Removing the power supply







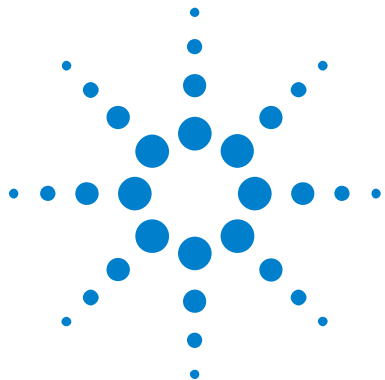

## 6 Replaceable Parts

This chapter describes how to order replaceable assemblies and parts for the Agilent 4000 X-Series oscilloscopes.

Diagrams and parts lists are included for assemblies and hardware that you can order.

Before working on the oscilloscope, read the safety summary at the back of this book.



## Ordering Replaceable Parts

### Listed Parts

To order a part in the parts list, quote the Agilent Technologies part number, indicate the quantity desired, and address the order to the nearest Agilent Technologies Sales Office. To find your nearest sales office go to **[www.agilent.com](http://www.agilent.com)**.

### Unlisted Parts

To order a part not listed in the parts list, include the instrument part number, instrument serial number, a description of the part (including its function), and the number of parts required. Address the order to the nearest Agilent Technologies Sales Office.

### Direct Mail Order System

Within the USA, Agilent Technologies can supply parts through a direct mail order system. There are several advantages to this system:

- Direct ordering and shipping from the Agilent Technologies parts center in California, USA.
- No maximum or minimum on any mail order. (There is a minimum amount for parts ordered through a local Agilent Technologies Sales Office when the orders require billing and invoicing.)
- Prepaid transportation. (There is a small handling charge for each order.)
- No invoices.

In order for Agilent Technologies to provide these advantages, please send a check or money order with each order.

Mail order forms and specific ordering information are available through your local Agilent Technologies Sales Office. Addresses and telephone numbers are located in a separate document shipped with the manuals.


## Exchange Assemblies

Some parts used in this instrument have been set up for an exchange program. This program allows the customer to exchange a faulty assembly with one that has been repaired, calibrated, and performance-verified by the factory. The cost is significantly less than that of a new part. The exchange parts have a part number in the form XXXXX-695XX.


After receiving the repaired exchange part from Agilent Technologies, a United States customer has 30 days to return the faulty assembly. For orders not originating in the United States, contact the local Agilent Technologies service organization. If the faulty assembly is not returned within 30 days, the customer will be charged an additional amount. The additional amount will be the difference in price between a new assembly and that of an exchange assembly.

## **Exploded Views**

The following exploded views provide a graphical representation of the oscilloscope at the time this manual was released. Not all parts are shown. Your parts may be slightly different than those shown. These views provide reference designator numbers that map to those used in the parts list table in this chapter.

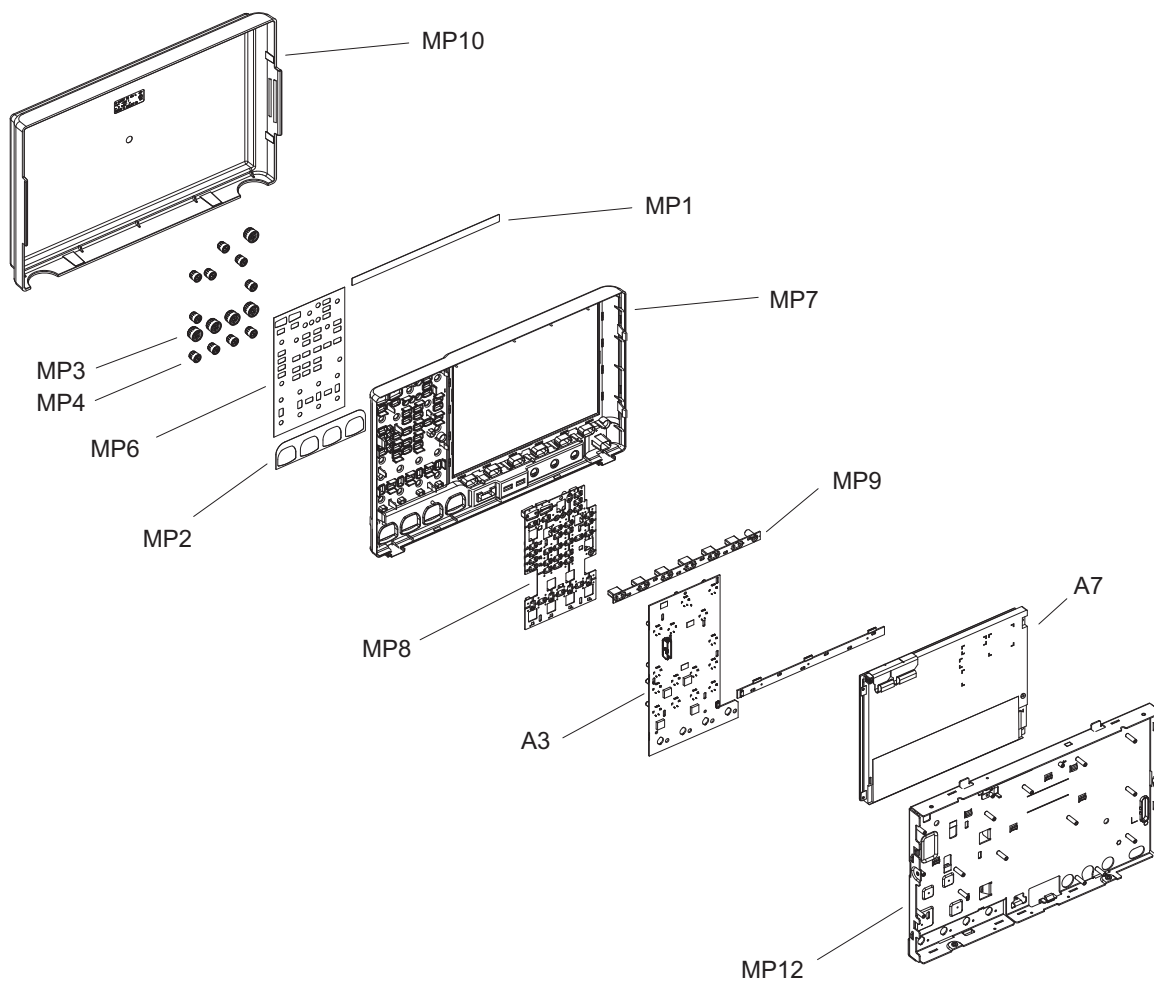

**Figure 53** Exploded View 1 of 2

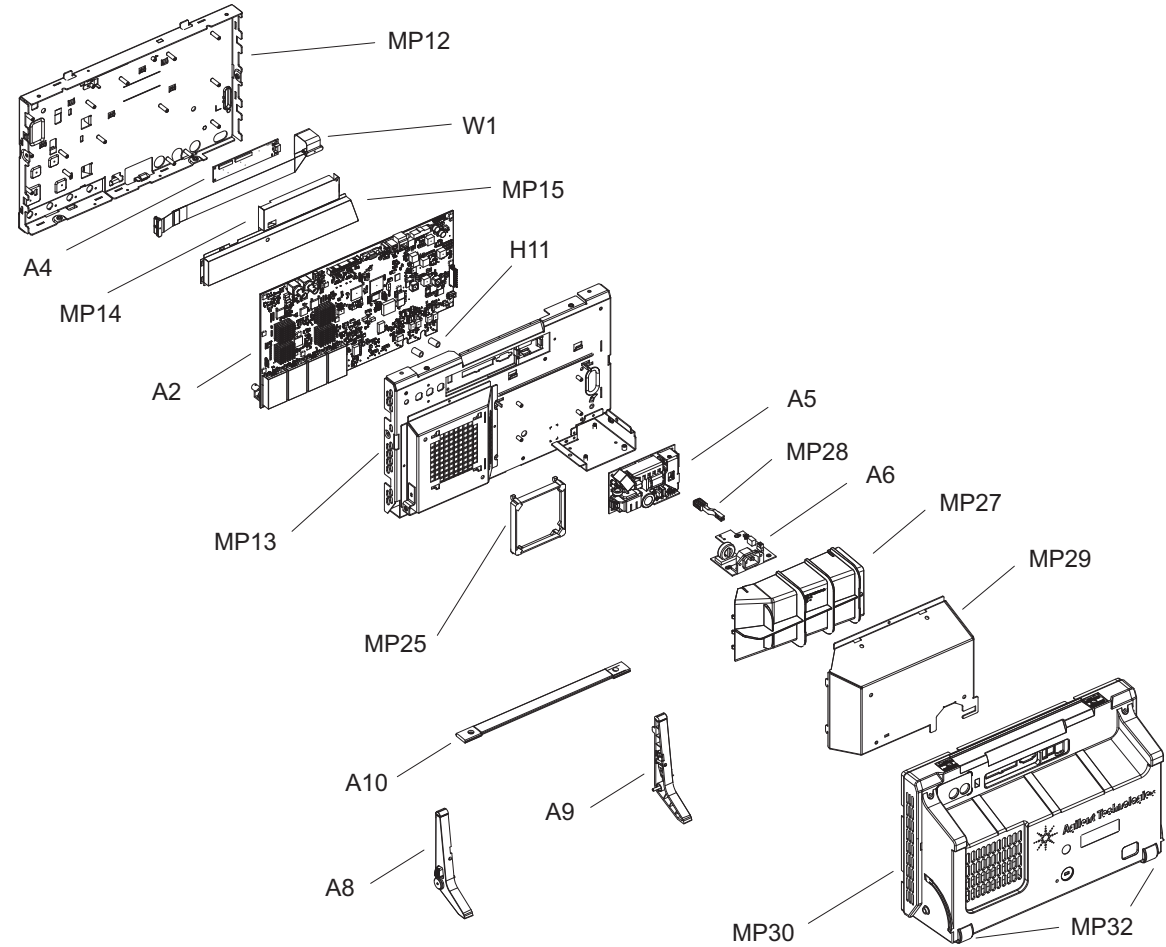

**Figure 54** Exploded View 2 of 2



## Replaceable Parts List

The information given for each part consists of the following:

- Reference designation.
- Agilent Technologies part number.
- Total quantity (Qty) in the instrument or on assembly.
- Description of the part.

**Table 19** Replaceable Parts

Ref Des	Agilent Part Number	Qty	Description
A2	54702-66601	1	Acquisition board 4000 X-Series 2-Ch - 200 MHz
A2	54702-69501	1	Acquisition board 4000 X-Series 2-Ch - 200 MHz (exchange)
A2	54702-69501-DEF	1	Acquisition board 4000 X-Series 2-Ch - 200 MHz (return)
A2	54703-66601	1	Acquisition board 4000 X-Series 4-Ch - 200 MHz
A2	54703-69501	1	Acquisition board 4000 X-Series 4-Ch - 200 MHz (exchange)
A2	54703-69501-DEF	1	Acquisition board 4000 X-Series 4-Ch - 200 MHz (return)
A2	54704-66601	1	Acquisition board 4000 X-Series 2-Ch - 350, 500 MHz
A2	54704-69501	1	Acquisition board 4000 X-Series 2-Ch - 350, 500 MHz (exchange)
A2	54704-69501-DEF	1	Acquisition board 4000 X-Series 2-Ch - 350, 500 MHz (return)
A2	54705-66601	1	Acquisition board 4000 X-Series 4-Ch - 350, 500 MHz
A2	54705-69501	1	Acquisition board 4000 X-Series 4-Ch - 350, 500 MHz (exchange)
A2	54705-69501-DEF	1	Acquisition board 4000 X-Series 4-Ch - 350, 500 MHz (return)
A2	54707-66601	1	Acquisition board 4000 X-Series 4-Ch - 1 GHz
A2	54707-69501	1	Acquisition board 4000 X-Series 4-Ch - 1 GHz (exchange)
A2	54707-69501-DEF	1	Acquisition board 4000 X-Series 4-Ch - 1 GHz (return)
A2	54709-66601	1	Acquisition board 4000 X-Series 4-Ch - 1.5 GHz
A2	54709-69501	1	Acquisition board 4000 X-Series 4-Ch - 1.5 GHz (exchange)

**Table 19** Replaceable Parts (continued)

Ref Des	Agilent Part Number	Qty	Description
A2	54709-69501-DEF	1	Acquisition board 4000 X-Series 4-Ch - 1.5 GHz (return)
A3	54709-66421	1	PCA - 2 Channel Keyboard 4000 X-Series
A3	54709-66422	1	PCA - 4 Channel Keyboard 4000 X-Series
A4	54709-66425	1	Touch Panel Controller board
A5	0950-5419	1	Power Supply Switching AC/DC 150W
A6	54709-66424	1	PCA - Line Filter Board
A7	54709-44701	1	Display assembly.
A8	54695-68708	1	Tilting Leg Right Assembly
A9	54695-68709	1	Tilting Leg Left Assembly
A10	54695-68711	1	Strap Handle Assembly
H2	54684-42202	2	Hole plug front panel (2 Channel only, not shown)
H5	2950-0054	2	Nut, 1/2-28 THD (secures BNCs, not shown)
H11	2190-0068	2	Washer-LK Intl T 1/2 In. .505 IN-ID (for BNCs, not shown)
H3	0515-0372	29	Screw-Machine w/Crest-Cup-Con-Wshr Pan-HD TORX-T10 M3X0.5 8mm-LG SST-300 Passivate (not shown)
H4	0515-0380	3	Screw-Machine W/Crest-Cup-Con-Washer Pan-HD Torx-T20 M4X0.7 10mm-LG SST-300 Passivated (not shown)
H5	0515-0658	2, 4	Screw-Machine w/Crest-Cup-Con-Wshr Pan-HD TORX-T6 M2X0.4 6mm-LG SST-300 Passivated (not shown, qty=2 for 2 Ch, qty=4 for 4 Ch)
H6	0515-0710	2	Screw-Machine 90-DEG-FLT-HD Torx-T20 M5X0.8 18mm-LG SST-300 Passivated (not shown)
H7	0515-2049	2	Screw-Machine w/Patch-LK 90-DEG-FLT-HD TORX-T20 M5X0.8 16mm-LG SST-300 Passivate (not shown)
H8	0515-2143	18	Screw-Machine w/Patch-LK Pan-HD TORX-T20 M4X0.7 6mm-LG SST-300 Passivate (not shown)
H9	3050-2320	2	Washer-Flat Non-metalic 10.5mm ID 20mm-OD 2mm-tHK Polyamide (not shown)

**Table 19** Replaceable Parts (continued)

Ref Des	Agilent Part Number	Qty	Description
H10	54695-26302	2	Screw-Shoulder (not shown)
H11	54709-24701	2	Spacer (at bottom of acquisition board)
MP1	54702-94301	1	ID Label DSO-X 4022A
MP1	54702-94302	1	ID Label MSO-X 4022A
MP1	54703-94301	1	ID Label DSO-X 4024A
MP1	54703-94302	1	ID Label MSO-X 4024A
MP1	54704-94301	1	ID Label DSO-X 4032A
MP1	54704-94302	1	ID Label MSO-X 4032A
MP1	54705-94301	1	ID Label DSO-X 4034A
MP1	54705-94302	1	ID Label MSO-X 4034A
MP1	54704-94303	1	ID Label DSO-X 4052A
MP1	54704-94304	1	ID Label MSO-X 4052A
MP1	54705-94303	1	ID Label DSO-X 4054A
MP1	54705-94304	1	ID Label MSO-X 4054A
MP1	54707-94301	1	ID Label DSO-X 4104A
MP1	54707-94302	1	ID Label MSO-X 4104A
MP1	54709-94301	1	ID Label DSO-X 4154A
MP1	54709-94302	1	ID Label MSO-X 4154A
MP2	54707-94304	1	BNC Label 2 Channel
MP2	54709-94304	1	BNC Label 4 Channel
MP3	75019-47401	5	Large Rotary Knob
MP4	75019-47402	10	Small Rotary Knob

## 6 Replaceable Parts

**Table 19** Replaceable Parts (continued)

Ref Des	Agilent Part Number	Qty	Description
MP6	54707-94303	1	Front Keyboard Label 2 Channel
MP6	54709-94303	1	Front Keyboard Label 4 Channel
MP7	54709-60201	1	Panel Bezel Assembly
MP8	54709-41901	1	Keyboard Keypad
MP9	54709-41903	1	Softkey Keypad
MP10	54709-44101	1	Front Cover
MP12	54709-00101	1	Front Deck
MP13	54709-00102	1	Rear Deck
MP14	54709-00602	1	Touch Screen Controller Shield
MP15	54709-00603	1	Keyboard Cable Shield
MP16	5067-4810	1	RFI Clip On Gasket (not shown)
MP17	54709-94327	1	Rear Label (not shown)
MP25	75019-44701	1	Fan Mount
MP26	3160-4365	1	Fan (not shown)
MP27	54709-47301	1	Air Duct
MP28	75019-43901	1	Extender, Power Switch
MP29	54695-00103	1	Power Supply Cover

**Table 19** Replaceable Parts (continued)

Ref Des	Agilent Part Number	Qty	Description
MP30	54709-64401	1	Cabinet Assembly
MP31	5042-6805	2	Cap - Strap Handle
MP32	54695-44701	2	Foot Pad
MP33	54709-94310	*	Overlay - Traditional Chinese, 4 channel (not shown)
MP33	54709-94309	*	Overlay - Traditional Chinese, 2 channel (not shown)
MP33	54709-94322	*	Overlay - Japanese, 4 channel (not shown)
MP33	54709-94321	*	Overlay - Japanese, 2 channel (not shown)
MP33	54709-94326	*	Overlay - Russian, 4 channel (not shown)
MP33	54709-94325	*	Overlay - Russian, 2 channel (not shown)
MP33	54709-94316	*	Overlay - French, 4 channel (not shown)
MP33	54709-94315	*	Overlay - French, 2 channel (not shown)
MP33	54709-94320	*	Overlay - Spanish, 4 channel (not shown)
MP33	54709-94319	*	Overlay - Spanish, 2 channel (not shown)
MP33	54709-94314	*	Overlay - German, 4 channel (not shown)
MP33	54709-94313	*	Overlay - German, 2 channel (not shown)
MP33	54709-94324	*	Overlay - Portuguese, 4 channel (not shown)
MP33	54709-94323	*	Overlay - Portuguese, 2 channel (not shown)
MP33	54709-94308	*	Overlay - Simplified Chinese, 4 channel (not shown)
MP33	54709-94306	*	Overlay - Simplified Chinese, 2 channel (not shown)
MP33	54709-94312	*	Overlay - Korean, 4 channel (not shown)
MP33	54709-94311	*	Overlay - Korean, 2 channel (not shown)
MP33	54709-94318	*	Overlay - Italian, 4 channel (not shown)
MP33	54709-94317	*	Overlay - Italian, 2 channel (not shown)

## 6 Replaceable Parts

**Table 19** Replaceable Parts (continued)

Ref Des	Agilent Part Number	Qty	Description
	54695-61605	1	Cable - Ground (not shown)
W1	54709-61601	1	Keyboard Cable
	54709-61603	1	Cable - DC (not shown)
	54709-61605	1	Cable - AC (not shown)
	54709-61606	1	Cable - Display (not shown)
	54709-61607	1	Cable - Backlight and Touch Screen Controller (not shown)
	54709-61608	1	Cable - Line Trigger (not shown)
	54709-61609	1	Cable - Remote (not shown)
	54709-61610	1	Soft Keyboard Cable (not shown)
Not shown	Power cord	0-1	Part number varies by country. Contact your local Agilent sales office for replacement.
Not shown	N2894A	*	Passive Probe 10:1, 700 MHz
Not shown	N2763-60001	*	Rack Mount Kit for 4000 X-Series Oscilloscope

\*Optional item.



## 7 Safety Notices

This apparatus has been designed and tested in accordance with IEC Publication 1010, Safety Requirements for Measuring Apparatus, and has been supplied in a safe condition. This is a Safety Class I instrument (provided with terminal for protective earthing). Before applying power, verify that the correct safety precautions are taken (see the following warnings). In addition, note the external markings on the instrument that are described under “Safety Symbols.”

### Warnings

Before turning on the instrument, you must connect the protective earth terminal of the instrument to the protective conductor of the (mains) power cord. The mains plug shall only be inserted in a socket outlet provided with a protective earth contact. You must not negate the protective action by using an extension cord (power cable) without a protective conductor (grounding). Grounding one conductor of a two-conductor outlet is not sufficient protection.

Only fuses with the required rated current, voltage, and specified type (normal blow, time delay, etc.) should be used. Do not use repaired fuses or short-circuited fuseholders. To do so could cause a shock or fire hazard.

If you energize this instrument by an auto transformer (for voltage reduction or mains isolation), the common terminal must be connected to the earth terminal of the power source.

Whenever it is likely that the ground protection is impaired, you must make the instrument inoperative and secure it against any unintended operation.

Service instructions are for trained service personnel. To avoid dangerous electric shock, do not perform any service unless qualified to do so. Do not attempt internal service or adjustment unless another person, capable of rendering first aid and resuscitation, is present.

Do not install substitute parts or perform any unauthorized modification to the instrument.

Capacitors inside the instrument may retain a charge even if the instrument is disconnected from its source of supply.

Do not operate the instrument in the presence of flammable gasses or fumes. Operation of any electrical instrument in such an environment constitutes a definite safety hazard.

Do not use the instrument in a manner not specified by the manufacturer.

## To clean the instrument

If the instrument requires cleaning: (1) Remove power from the instrument. (2) Clean the external surfaces of the instrument with a soft cloth dampened with a mixture of mild detergent and water. (3) Make sure that the instrument is completely dry before reconnecting it to a power source.



## Safety Symbols



Instruction manual symbol: the product is marked with this symbol when it is necessary for you to refer to the instruction manual in order to protect against damage to the product.



Hazardous voltage symbol.



Earth terminal symbol: Used to indicate a circuit common connected to grounded chassis.



# Index

## A

accuracy, DC vertical gain, [33](#)  
accuracy, dual cursor, [39](#)  
Agilent  
    contact information, [19](#)  
assemblies  
    exchange, [137](#)  
    replacing, [99](#)

## B

bandwidth  
    verification, [45](#)

## C

cabinet, removing, [102](#), [115](#)  
calibration, [67](#)  
cautions, [83](#), [86](#), [89](#), [100](#), [101](#), [132](#)  
certificate of calibration, [67](#)  
channels  
    problem solving, [74](#)  
cleaning the instrument, [148](#)  
compensating your probe, [80](#)  
connecting for threshold test, [30](#)  
contact Agilent, [19](#)

## D

DC vertical gain accuracy, [33](#)  
digital channel  
    testing, [27](#)  
    verification, [28](#)  
digital channels, [4](#)  
display  
    assembly removal, [116](#)  
    problem solving, [72](#)  
dual cursor accuracy, [39](#)

## E

equipment, test, [23](#)  
exchange assemblies, [137](#)  
exploded view  
    front panel, [139](#)  
    instrument, [140](#)

## F

fan assembly removal, [125](#)  
front panel  
    exploded view, [139](#)  
    removal, [105](#)  
front panel self test, [76](#)

## H

handle, removing, [103](#)  
hardware self test, [76](#)

## I

internal self-tests, [76](#)

## K

keyboard, removal, [109](#)

## L

LED, startup sequence, [75](#)  
lid, removal, [104](#)  
line filter board removal, [129](#)

## M

making test connector, [25](#)  
measurement category, [15](#)  
    definitions, [15](#)

## N

NIST, [67](#)

## O

ordering parts, [136](#)  
overvoltage category, [17](#)

## P

parts  
    list, [141](#)  
    ordering, [136](#)  
passive probe, [146](#)  
perform user cal, [67](#)  
performance  
    test record, [61](#)  
pollution degree, [17](#)  
pollution degree, definitions, [17](#)  
power  
    cord list, [146](#)  
    supply removal, [132](#)  
    supply shield, removal, [126](#)  
    test points, [87](#), [88](#)  
power requirements, [14](#)  
probe  
    N2894A, [146](#)  
probes  
    compensating, [80](#)  
    problem solving, [71](#)

## R

record tests, [61](#)

## Index

remove, 129  
    cabinet, 102, 115  
    display assembly, 116  
    fan assembly, 125  
    front panel, 105  
    handle, 103  
    keyboard, 109  
    power supply, 132  
    power supply shield, 126  
    storage lid, 104  
replaceable parts list, 141  
replacing assemblies, 99  
run self-tests, 76

## S

safety notice, 65, 83  
self test, front panel, 76  
self test, hardware, 76  
self-alignment, user cal, 67  
self-tests, 76  
specifications, 18  
startup sequence, 75  
status, User Cal, 70  
status, user cal, 67  
storage lid, removal, 104

## T

test  
    connector, constructing, 25  
    digital channels, 27  
    equipment, 23  
    points, power, 87, 88  
    record, 61  
threshold  
    accuracy, digital channel, 28  
    test diagram, 30  
trace display  
    problem solving, 73  
transient withstand capability, 16  
trigger  
    sensitivity, 53  
    sensitivity, external, 57  
    sensitivity, internal, 54

troubleshooting, 71  
    equipment required, 82  
    fan, 97  
    flowchart, 86  
    power supply, 89

## U

user calibration, 67

## V

verify  
    bandwidth, 45  
    digital channel, 28  
    trigger, 53

## W

warnings, 100, 127  
warranted specifications, 18